

---

# **AXONICA<sup>®</sup> 3.1**

## User Guide





Copyright (C) 2018, Team aXonica

**Team aXonica, 2017-2018**

Bilal Wajid  
Hasan Iqbal  
Momina Jamil

**aXonica 3.1 User Guide**

This program is free software: you can redistribute it and/or modify it under the terms of the **GNU General Public License** as published by the Free Software Foundation, either **version 3** of the License, or (at your option) any later version.

This program is distributed in the hope that it will be useful, but WITHOUT ANY WARRANTY; without even the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR PURPOSE. See the **GNU General Public License** for more details.

You should have received a copy of the GNU General Public License along with this program. If not, see:

<https://www.gnu.org/licenses/>

**Disclaimer of Warranty:**

THERE IS NO WARRANTY FOR THE PROGRAM, TO THE EXTENT PERMITTED BY APPLICABLE LAW. EXCEPT WHEN OTHERWISE STATED IN WRITING THE COPYRIGHT HOLDERS AND/OR OTHER PARTIES PROVIDE THE PROGRAM "AS IS" WITHOUT WARRANTY OF ANY KIND, EITHER EXPRESSED OR IMPLIED, INCLUDING, BUT NOT LIMITED TO, THE IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE. THE ENTIRE RISK AS TO THE QUALITY AND PERFORMANCE OF THE PROGRAM IS WITH YOU. SHOULD THE PROGRAM PROVE DEFECTIVE, YOU ASSUME THE COST OF ALL NECESSARY SERVICING, REPAIR OR CORRECTION.

---

**Limitation of Liability:**

IN NO EVENT UNLESS REQUIRED BY APPLICABLE LAW OR AGREED TO IN WRITING WILL ANY COPYRIGHT HOLDER, OR ANY OTHER PARTY WHO MODIFIES AND/OR CONVEYS THE PROGRAM AS PERMITTED ABOVE, BE LIABLE TO YOU FOR DAMAGES, INCLUDING ANY GENERAL, SPECIAL, INCIDENTAL OR CONSEQUENTIAL DAMAGES ARISING OUT OF THE USE OR INABILITY TO USE THE PROGRAM (INCLUDING BUT NOT LIMITED TO LOSS OF DATA OR DATA BEING RENDERED INACCURATE OR LOSSES SUSTAINED BY YOU OR THIRD PARTIES OR A FAILURE OF THE PROGRAM TO OPERATE WITH ANY OTHER PROGRAMS), EVEN IF SUCH HOLDER OR OTHER PARTY HAS BEEN ADVISED OF THE POSSIBILITY OF SUCH DAMAGES.



# Contents

About This Guide.....	1
Abstract.....	1
Purpose.....	1
Intended Audience .....	1
Required Background.....	1
Updates from Previous Release.....	1
How This Guide is Organized.....	2
Contact.....	2
aXonica Package.....	3
Installing aXonica software .....	3
Software Installation Process.....	5
Installing aXonica software .....	5
Microsoft Windows Support .....	13
Installing aXonica software on Microsoft Windows .....	13
aXonica-UI Application .....	23
In this Chapter.....	23
Introduction.....	23
Configuration of aXonica-UI Application.....	24
Data Acquisition and Software Pipeline.....	26
In this Chapter.....	26
Introduction.....	26
Image Data Acquisition from MRI .....	27
Software Image Pipeline .....	29
Tutorial Datasets .....	30
Downloading aXonica Tutorial Datasets .....	30
Image Preprocessing.....	34
In this Chapter.....	34
Introduction.....	34
Image Visualization.....	35
ImageVis3D [I].....	35

Image Enhancement.....	37
AMIDE [2].....	37
3D Image Reconstruction.....	38
ITK-SNAP [3].....	38
Image Processing.....	40
In this Chapter.....	40
Introduction.....	40
Image Segmentation.....	41
BrainSeg3D [4].....	41
Image Registration.....	43
Mango [5].....	43
Image Surfacing.....	45
BrainVISA [6].....	45
Structural Analysis.....	48
In this Chapter.....	48
Introduction.....	48
Brain Connectomics.....	49
medInria [7].....	50
3D Image Analysis.....	52
ImageJ [8].....	52
Morphometric Analysis.....	53
FracLac [9].....	53
Image Classification.....	56
LA-iMageS [10].....	56
Image Mapping.....	58
MRI Processor [11].....	58
Data Management and Annotation.....	61
In this Chapter.....	61
Introduction.....	61
Image Format Management.....	62
MRICron [12].....	62
DICOM File Management.....	63
WEASIS [13].....	63
Post Installation.....	65
Post Installation of aXonica software.....	65
Anatomist.....	65
Dipy.....	65

---

---

BoneJ.....	65
ImageVIS 3D.....	65
AMIDE.....	65
PID.....	65
TORTOISE.....	66
normalizeFOV.....	66
FastFilter3D.....	66
ITK-Snap.....	66
HeteroscedasticfMRI.....	66
pyClusterROI.....	66
BrainSeg3D .....	66
CMP-BIA.....	66
ITK.....	66
E-Snake.....	67
mincBEAST .....	67
Segmentator.....	67
Mango.....	67
bUnwarpJ .....	67
BrainVISA.....	67
The Virtual Brain .....	67
CBS Tools.....	67
HegaEtAI2017.....	68
NEURON.....	68
SliceMap.....	68
Nengo.....	68
neuroConstruct.....	68
MedInria.....	68
Brian.....	68
Time Domain decoding .....	68
Fraclac.....	68
Gwyddion.....	69
LA-iMageS.....	69
Brainstorm .....	69
ImageJ .....	69
Nipype.....	69
MRI Processor.....	69
NiBabel.....	69
MRIcron.....	69

---

MIPAV .....	69
WEASIS .....	70
Software Uninstallation Process .....	71
Uninstalling aXonica software .....	71
References.....	76

# About This Guide

---

## Abstract

aXonica is a free, easily distributable software installation package for freeware tools related to Bioimaging especially in the field of Magnetic Resonance Imaging (MRI).

aXonica is a Linux-based installation package that downloads and installs the tools related to Magnetic Resonance Imaging (MRI) as per user requirement. It also downloads the dependencies prior to the tools for which they are required. aXonica also provides an interactive user interface which is easily understandable for users unaware of UNIX-shell language.

---

## Purpose

The basic purpose of this document is to provide the complete description about the installation of aXonica shell file as well as the complete working of aXonica-UI. This document also explains the complete pipeline of MRI tools; also shows how this pipeline is implemented in the installation package.

---

## Intended Audience

This document is intended for users who wish to install the tools related to Magnetic Resonance Imaging on Linux-based operating system.

---

## Required Background

Team aXonica has made every attempt to make this a step by step guide. However, some familiarity with Linux operating system as well as software and hardware requirements of aXonica are assumed.

---

## Updates from Previous Release

Added Microsoft Windows Virtualization Support Chapter which enables the user to work on a Windows platform.

---

## How This Guide is Organized

This guide is organized into sections grouped according to the intended use by the user:

- About This Guide (Chapter 1) describes this document's purpose and intended audience.
  - aXonica Package (Chapter 2) contains a table which enlists all the included software.
  - Software Installation Process (Chapter 3) describes how to install aXonica.
  - Microsoft Windows Support (Chapter 4) describes a step by step procedure to configure a virtual OS for aXonica on Windows.
  - aXonica-UI (Chapter 5) describes the step by step working of user interface.
  - Data Acquisition and Software Pipeline (Chapter 6) describes the acquisition of raw data image of MRI machine and implementation of software pipeline of Biological image development.
  - Tutorial Datasets (Chapter 7) describes how to download tutorial datasets.
  - Image Preprocessing (Chapter 8) describes the step by step solution to recommended software.
  - Image Processing (Chapter 9) describes the step by step solution to recommended software.
  - Structural Analysis (Chapter 10) describes the step by step solution to recommended software.
  - Data Management & Annotation (Chapter 11) describes the step by step solution to recommended software.
  - Post Installation (Chapter 12) describes the working of each software included in the package.
  - Software Uninstallation Process (Chapter 13) describes how to install aXonica or any specific tools.
- 

## Contact

For any further queries and suggestions, contact us at: [hasaniqbal777@gmail.com](mailto:hasaniqbal777@gmail.com) or [mominaj05@gmail.com](mailto:mominaj05@gmail.com)

# aXonica Package

## Installing aXonica software

List of software tools and plugins included in aXonica package are shown in the following table. Software which are recommended are also mentioned in the table. Download size and version of each software is mentioned for user convenience:

Sr. no.	Software Name	Size(MB)	Software Version	Recommended
<b>Preprocessing</b>				
Image Visualization				
1	Anatomist	N-A	4.5.0	
2	Dipy	30	0.15.0	
3	BoneJ	0.7	1.0.0	
4	ImageVIS3D	30.9	3.1.0	✓
Image Enhancement				
5	AMIDE	40	1.0.4	✓
6	PID	5.1	1.12.0	
7	TORTOISE	751.5	3.1.0	
8	NormalizeFOV	732	1.2.0	
9	FastFilter3D	0.3	1.0.0	
3D Image Reconstruction				
10	ITK-SNAP	45.2	3.6.0	✓
<b>Processing</b>				
Image Segmentation				
11	Heteroscedastic	0.19	1.0.0	
12	pyClusterROI	0.4	1.0.0	
13	BrainSeg3D	18.5	1.0.0	✓
14	CMP-BIA	0.1	0.3.0	
15	ITK	88	2.4.7	
16	E-Snake	0.1	1.0.0	
17	mincBEAST	1.96	1.15.0	
18	Segmentator	1.7	0.19.1	

Image Registration				
19	Mango	67	1.7.0	✓
20	bunwarpJ	0.2	2.6.5	
Surfacing				
21	BrainVISA	1100	4.5.0	✓
Structural Analysis				
Brain Connectomics				
22	The Virtual Brain	670	1.5.4	
23	CBS Tools	147	3.1.0	
24	HagaEtAl	0.1	1.0.0	
25	NEURON	8	7.5.0	
26	SliceMap	44	1.0.0	
27	Nengo	38	2.8.0	
28	neuroConstruct	0.4	1.6.0	
29	MedInria	58	2.2.3	✓
30	BRIAN	68	2.0.2	
31	Timedomain decoding	0.5	1.0.0	
Morphometric Analysis				
32	FracLac	13.8	1.0.0	✓
33	Gwyddion	50	2.44.1	
Classification				
34	LA-iMages	298	1.1.5	✓
3D Image Analysis				
35	Brainstorm	76.6	3.0.0	
36	ImageJ	137.7	2.0.0	✓
37	nipype	48	1.1.7	
Image Mapping				
38	MRI-Processor	0.5	1.0.0	✓
Image Management and Annotation				
Format Management				
39	NiBabel	48	2.3.1	
40	MRIcron	24.5	1.0.0	✓
DICOM File Management				
41	MIPAV	134	8.0.2	
42	WEASIS	19	3.0.4	✓

# Software Installation Process

## Installing aXonica software

For downloading of aXonica, visit its website:

<https://github.com/hasaniqbal777/aXonica-bin>

1. aXonica can also be downloaded from the following command through git:  
`$ git clone https://github.com/hasaniqbal777/aXonica-bin`
2. Now run the following commands on terminal:  
`$ chmod +x aXonica_setup_enUS`  
`$ sudo bash aXonica_setup_enUS`
3. Installation wizard of aXonica will start.

**NOTE** *Installation of aXonica require a proper internet connection to proceed, otherwise the installation terminates.*

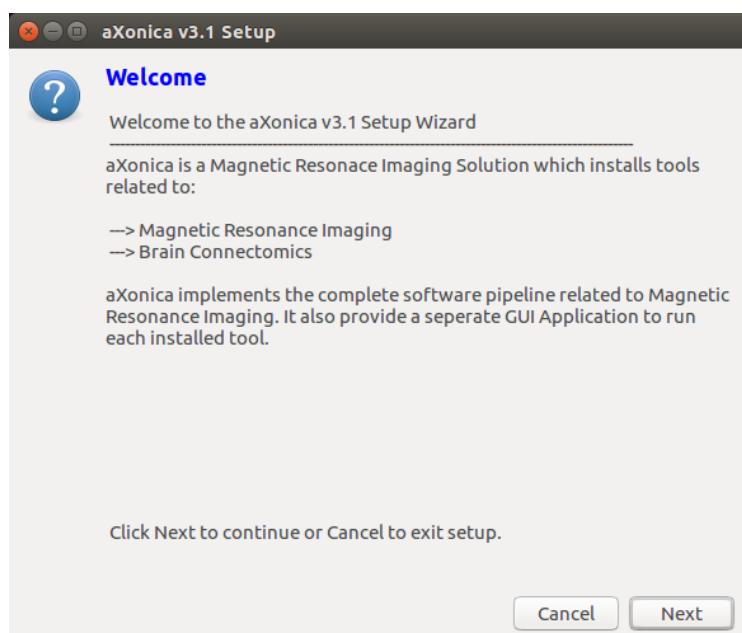


Figure 1: Installation Welcome screen

4. Click **Next** to proceed and confirm the Internet connection.

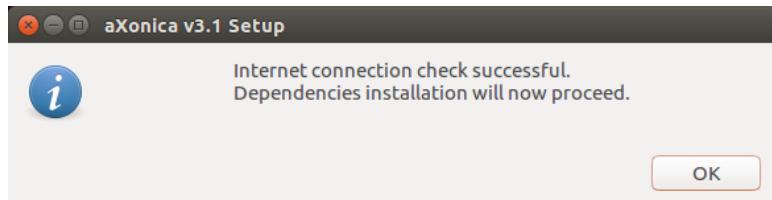


Figure 2: Internet check successful screen

5. Installation is **terminated** if there is **no internet**.

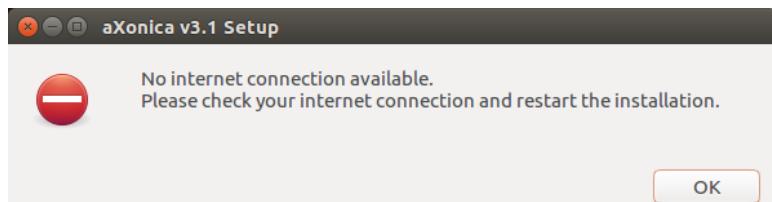


Figure 3: Internet check failed screen

6. **Dependencies** related to software will start installing.

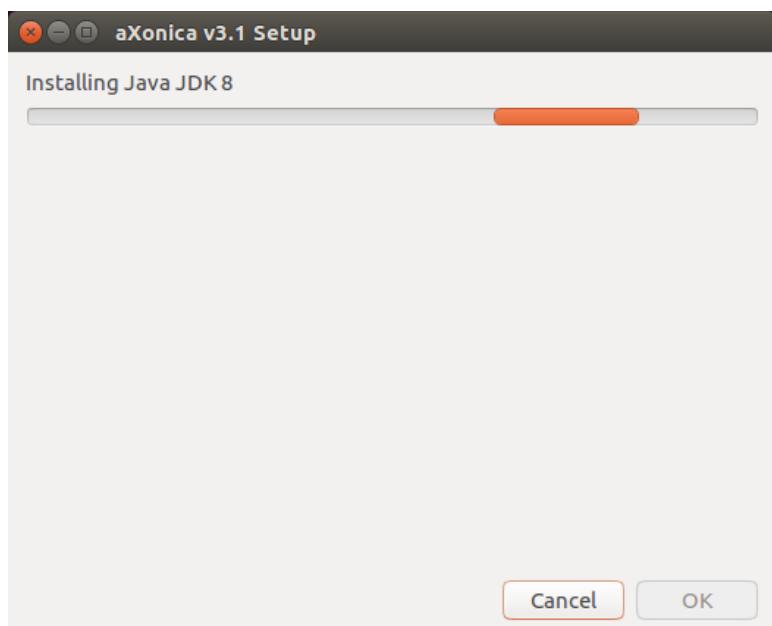


Figure 4: Dependencies Installation in progress screen

7. All the dependencies which are installed are displayed at the end of the installation. Click **Ok**.

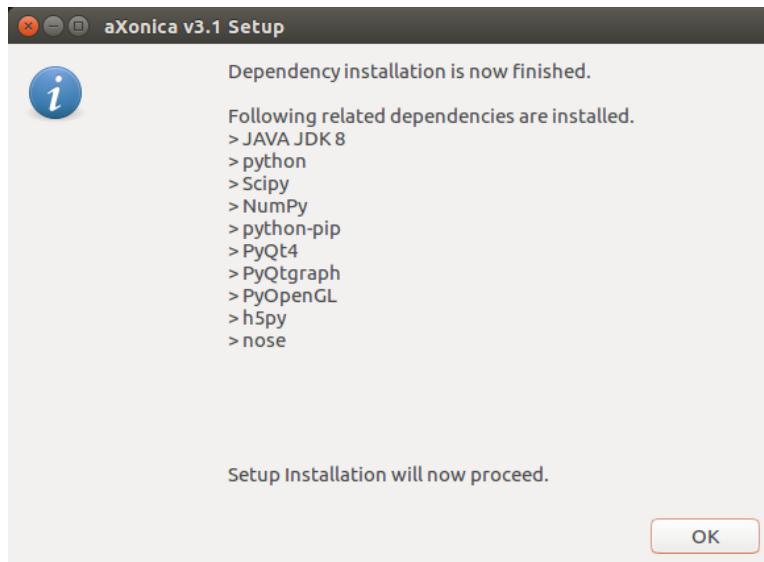


Figure 5: Dependencies Installation finished screen

8. aXonica pipeline detail screen is **displayed**. It has the information about which software you want to install. Click **Next** to Proceed.

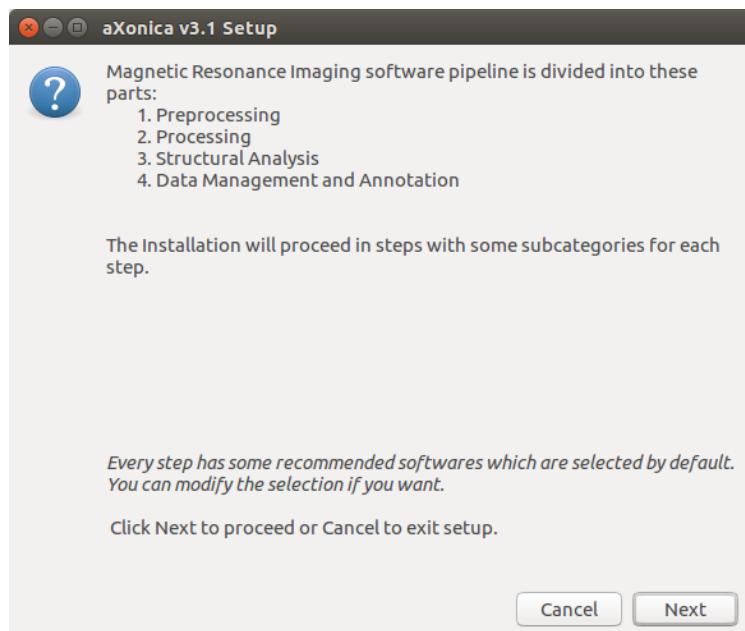


Figure 6: Pipeline information screen

9. Selection screen for Preprocessing tools is displayed. Select the required tools and Click **Next**.

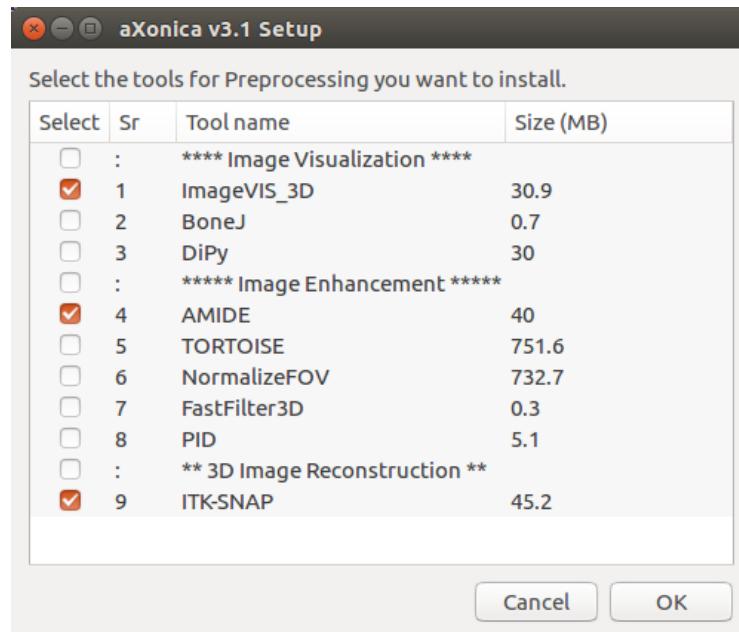


Figure 7: Preprocessing tools selection screen

10. Selection screen for Processing tools is displayed. Select the required tools and Click **Next**.

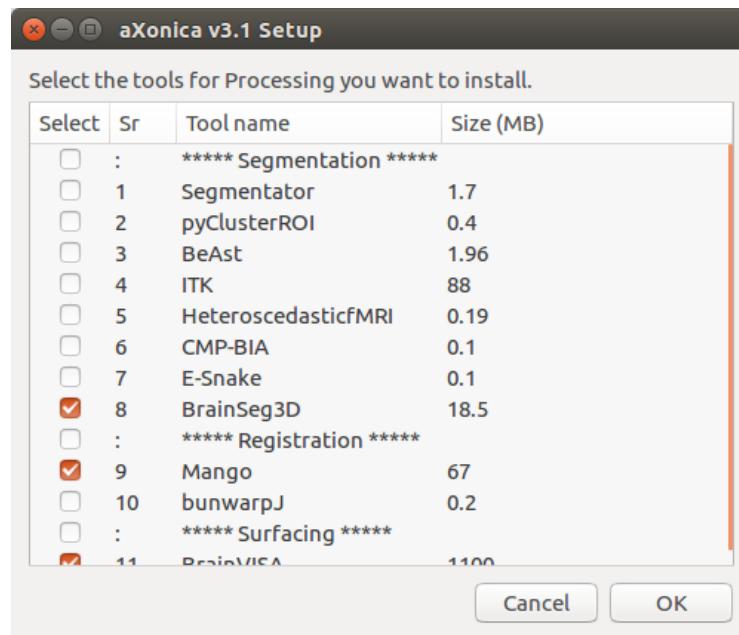


Figure 8: Processing tools selection screen

11. Selection screen for Structural Analysis tools is displayed. Select the required tools and Click **Next**.

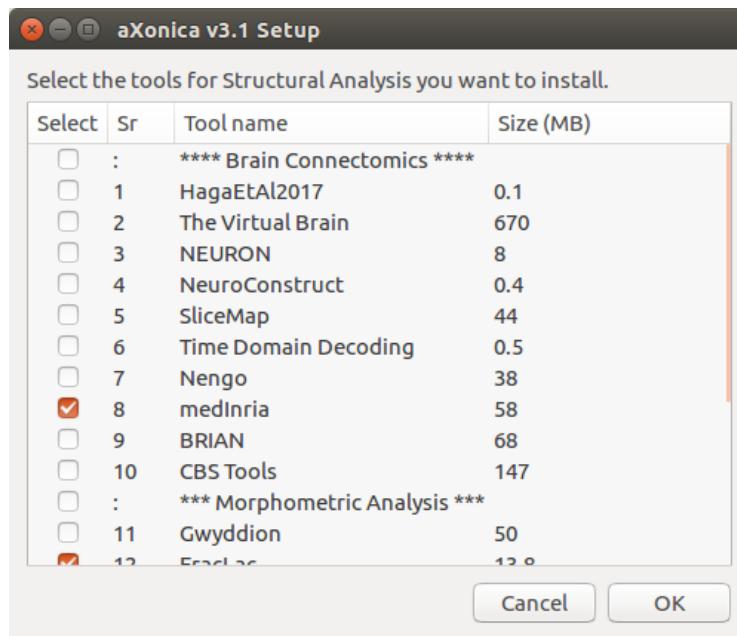


Figure 9: Structural Analysis tools selection screen

12. Selection screen for Data Management and Annotation tools is displayed. Select the required tools and Click **Next**.

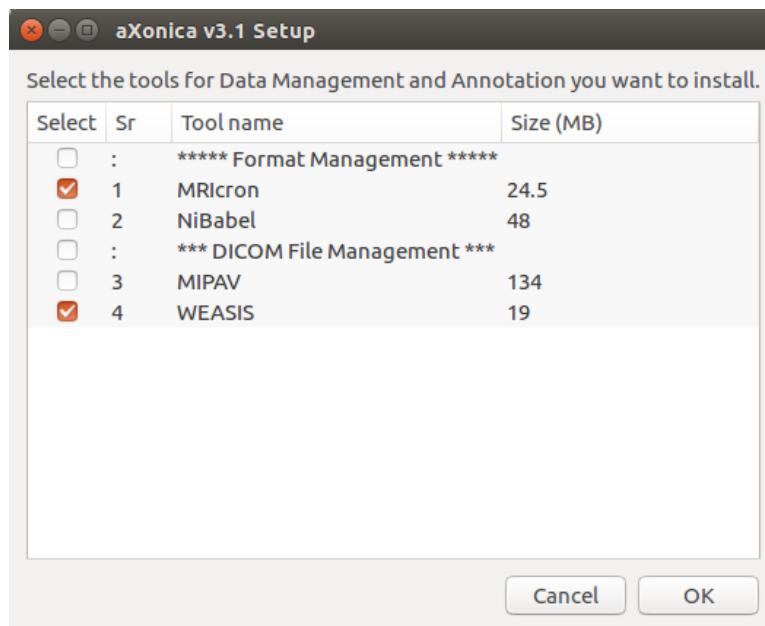


Figure 10: Data Management and Annotation selection screen

13. Click **Next** to proceed with the installation.

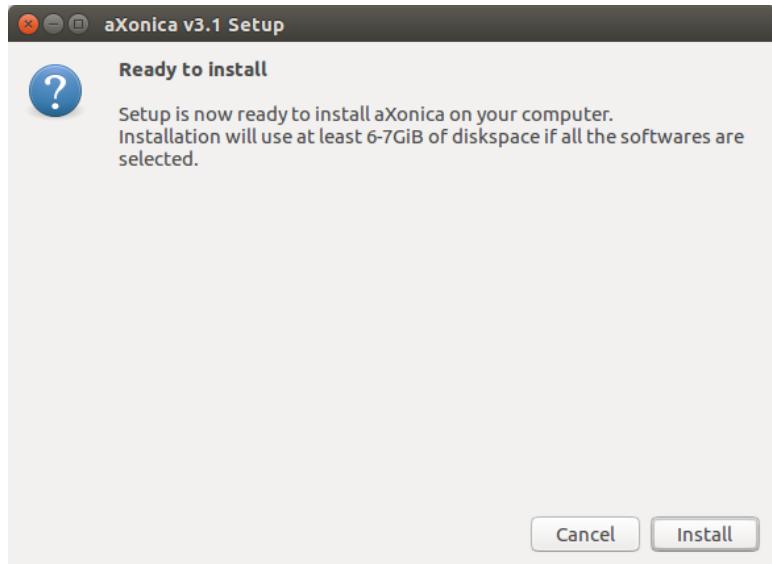


Figure 11: Ready to Install screen

14. Tools installation will **continue**. Some tools **install** as a standalone installation and will be called automatically.


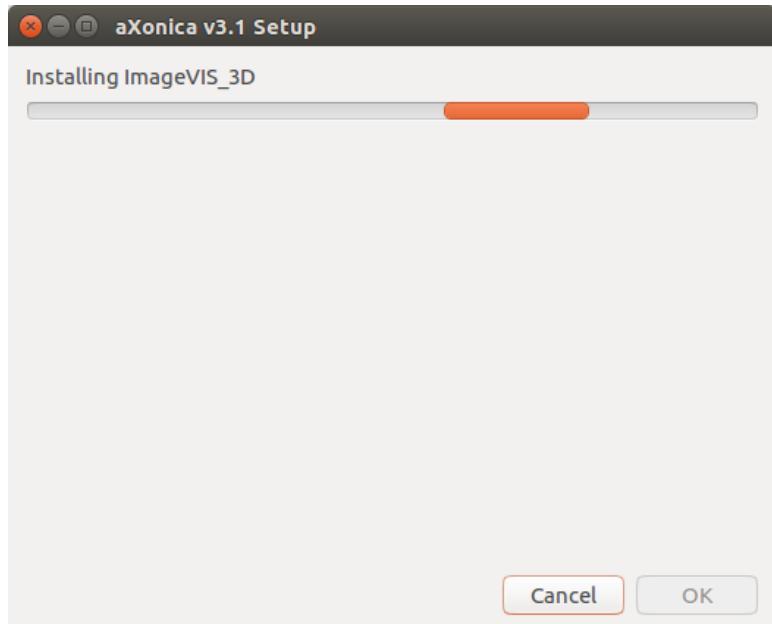


Figure 12: Tools installation progress screen

15. Software will be **download** from the aXonica repository online. Downloading time depend upon your Internet.

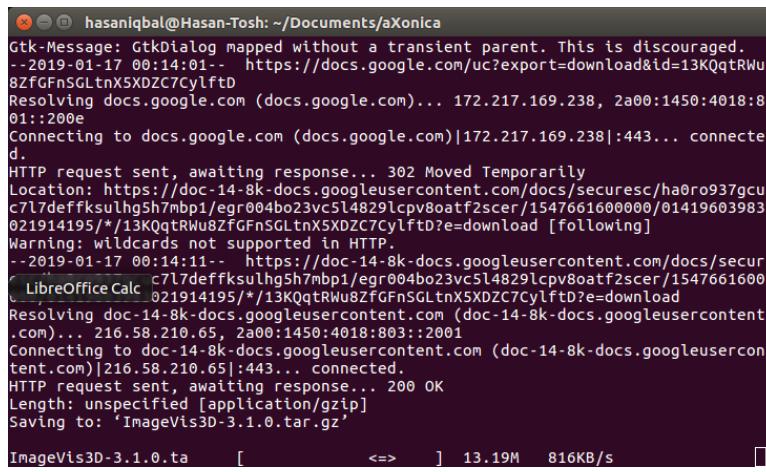

```
hasaniqb@Hasan-Tosh:~/Documents/aXonica
Gtk-Message: GtkDialog mapped without a transient parent. This is discouraged.
--2019-01-17 00:14:01-- https://docs.google.com/uc?export=download&id=13KQqtRWu
8ZfGFnSGLtnX5XDZC7CylfD
Resolving docs.google.com (docs.google.com)... 172.217.169.238, 2a00:1450:4018:8
01::200e
Connecting to docs.google.com (docs.google.com)|172.217.169.238|:443... connecte
d.
HTTP request sent, awaiting response... 302 Moved Temporarily
Location: https://doc-14-8k-docs.googleusercontent.com/docs/securesc/ha0ro937gcu
c7l7deffksulg5h7mbp1/egr004bo23vc5l4829lcpv8oatf2scer/1547661600000/01419603983
021914195/*/13KQqtRWu8ZfGFnSGLtnX5XDZC7CylfD?e=download [following]
Warning: wildcards not supported in HTTP.
--2019-01-17 00:14:11-- https://doc-14-8k-docs.googleusercontent.com/docs/secur
[ LibreOffice Calc c7l7deffksulg5h7mbp1/egr004bo23vc5l4829lcpv8oatf2scer/1547661600
021914195/*/13KQqtRWu8ZfGFnSGLtnX5XDZC7CylfD?e=download
Resolving doc-14-8k-docs.googleusercontent.com (doc-14-8k-docs.googleusercontent
.com)... 216.58.210.65, 2a00:1450:4018:803::2001
Connecting to doc-14-8k-docs.googleusercontent.com (doc-14-8k-docs.googleusercontent
.com)|216.58.210.65|:443... connected.
HTTP request sent, awaiting response... 200 OK
Length: unspecified [application/gzip]
Saving to: 'ImageVls3D-3.1.0.tar.gz'

ImageVls3D-3.1.0.ta [=====] 13.19M 816KB/s

```

Figure 13: Tools downloading in background

16. Setup will now install aXonica\_UI Application on your computer. Click **Next** to proceed.

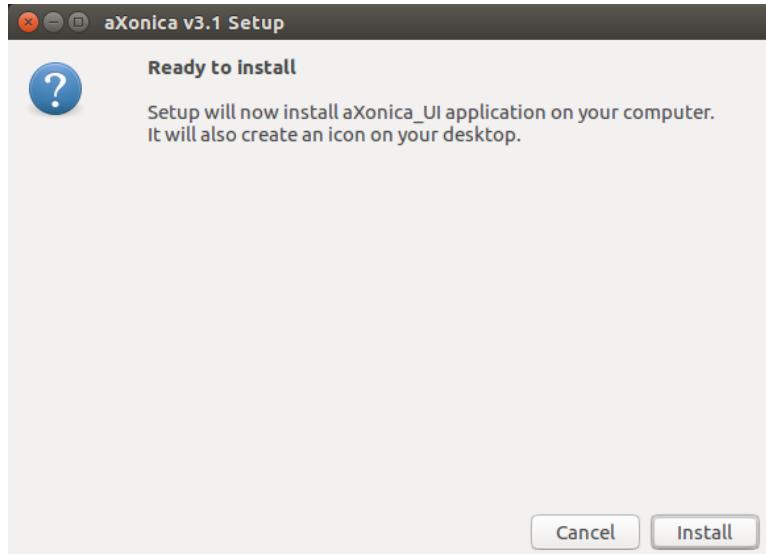


Figure 14: Ready to Install aXonica\_UI screen

17. Installation of aXonica\_UI will now proceed.

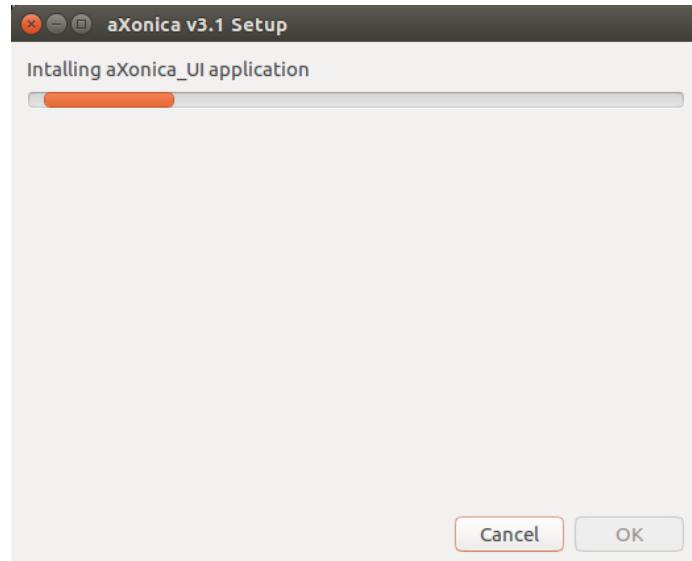


Figure 15: aXonica\_UI Installation progress screen

18. Installation of aXonica is now finished. Click **Finish** to use the Tools.

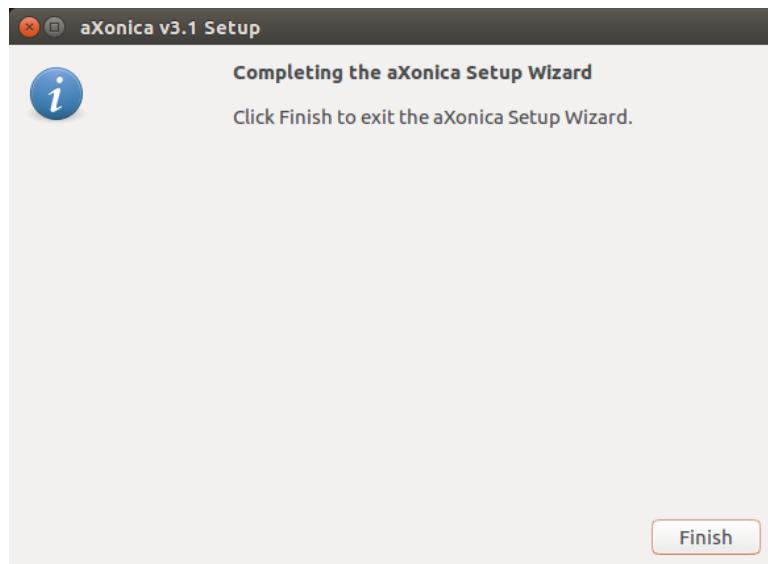


Figure 16: Finalizing Installation screen

# Microsoft Windows Support

## Installing aXonica software on Microsoft Windows

aXonica is software package for Linux, specifically Ubuntu users, but Windows users can use aXonica using a Virtual Machine. Follow these steps to initialize a Virtual Machine on your Windows host.

**NOTE** aXonica support a 64bit Windows host. Microsoft Windows XP and Vista support is discontinued.

1. First you have to **download** the image file of the required OS (Ubuntu) from its website or use the following link:  
<http://repo.isra.edu.pk/ubuntu-release/18.04.2/ubuntu-18.04.2-desktop-amd64.iso>
2. Download the virtualization software (Oracle VM VirtualBox) from the following link:  
<https://download.virtualbox.org/virtualbox/6.0.10/VirtualBox-6.0.10-132072-Win.exe>
3. **Install** this software in Windows OS host.
4. Start Virtual Box, and click on the **New** symbol.

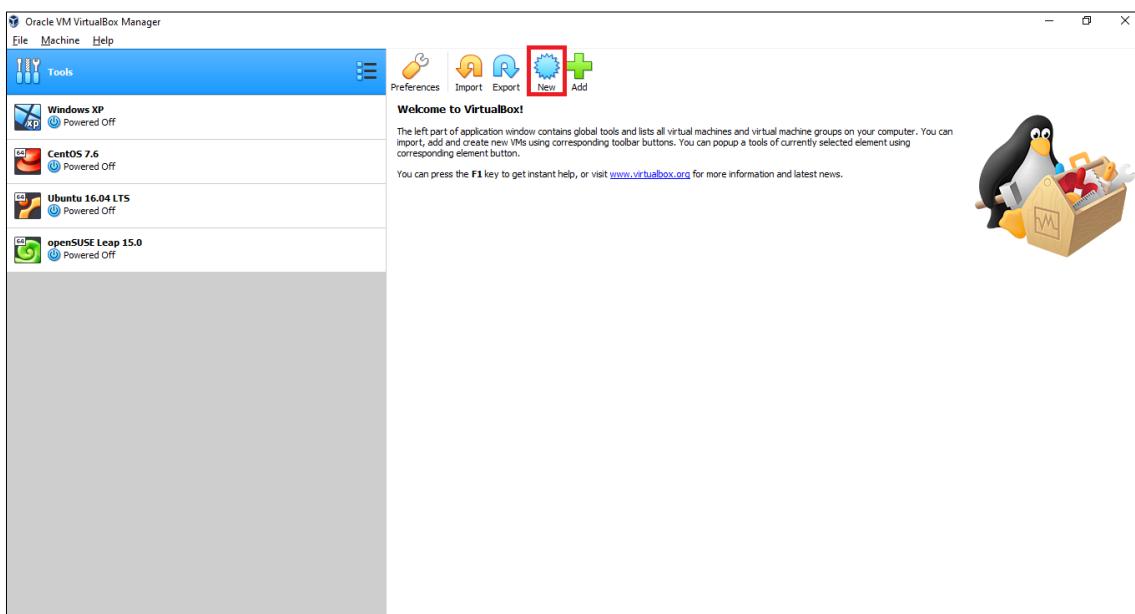


Figure 1: New Virtual OS

5. Give the virtual OS a relevant **Name**. Select the **Type** (Linux) and **Version** (Ubuntu 64-bit) and Click **Next**.

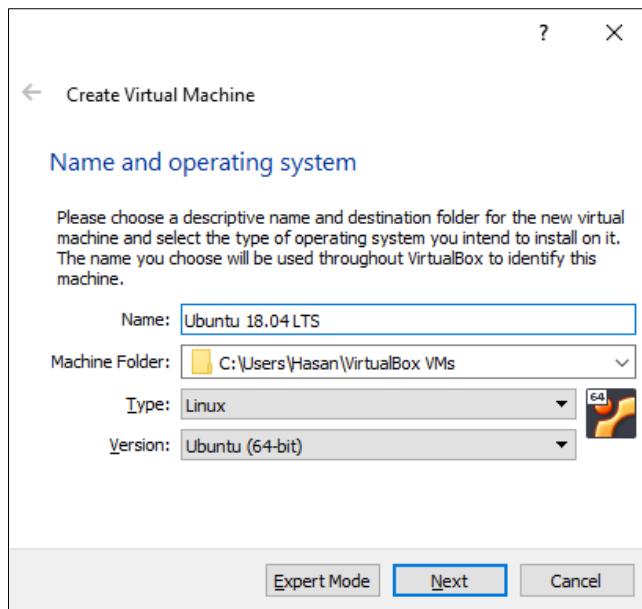


Figure 2: Assigning Name

6. **Allocate** RAM to the virtual OS. Following system has 8GB of RAM so 2GB of RAM is allocated. You can use more RAM if your system has enough extra RAM.

**NOTE** Allocate about half of the RAM to the virtual OS. Click **Next**.

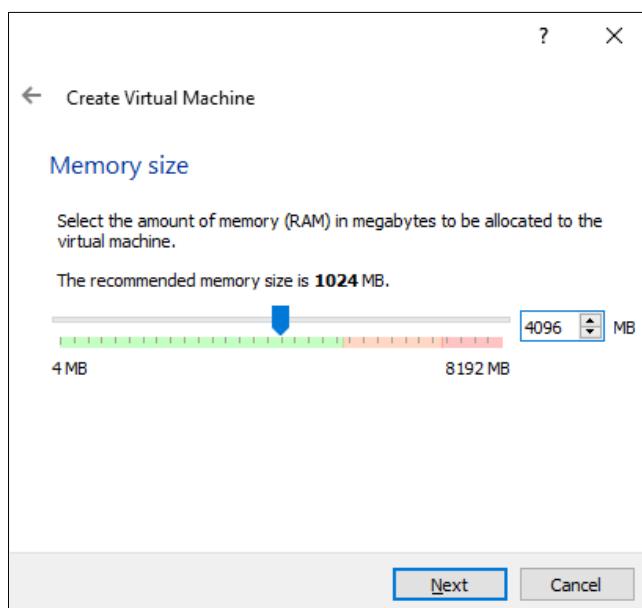


Figure 3: Memory Allocation

7. **Create a virtual disk.** This works as the hard disk of the virtual Linux system. This is where the virtual system will store its files. Click **Create**.

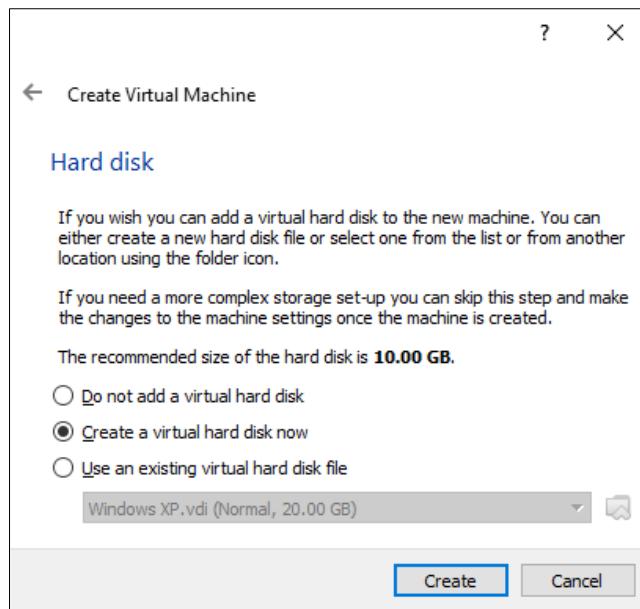


Figure 4: Creating a virtual hard disk

8. Select VDI file type here (recommended). Click **Next**.

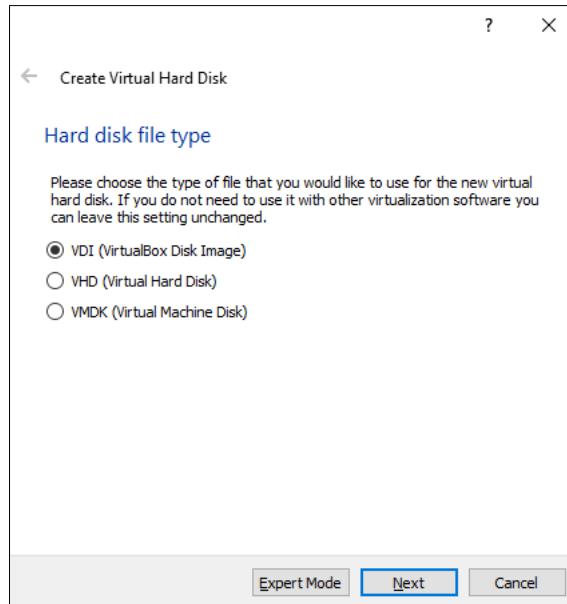


Figure 5: Hard disk file type

9. You can choose either of Dynamically allocated or Fixed size option for creating the virtual hard disk. Choose Dynamically allocated. (recommended). Click **Next**.

**NOTE** *Dynamic allocation is allocated as time passes and data is increased whereas fixed is allocated instantly.*

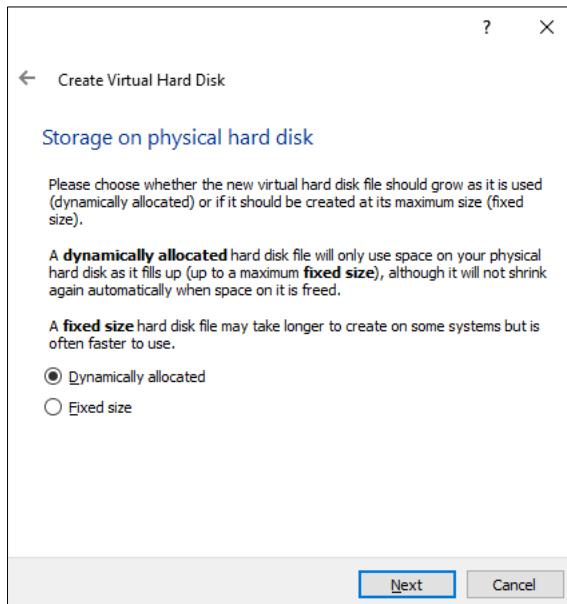


Figure 6: Storage type

10. Select VDI file type here (recommended) and Select the Hard Disk size. (recommended size: 100 GB). Click **Create**.

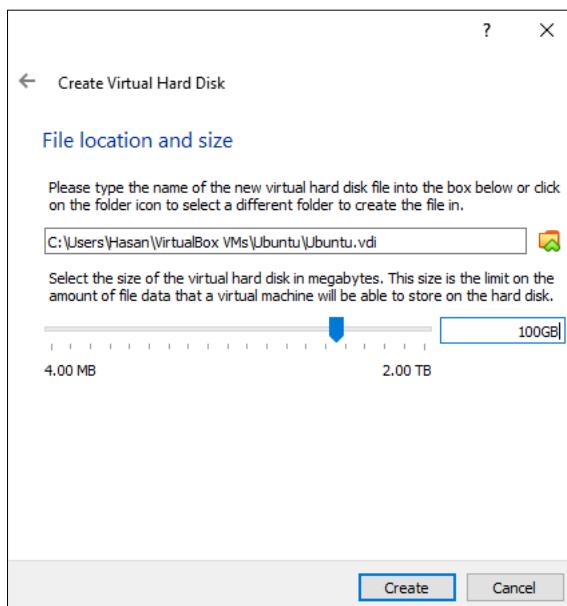


Figure 7: Hard disk size

11. Click **Next**. Now, Select Settings to assign the image file of respective OS to VB.

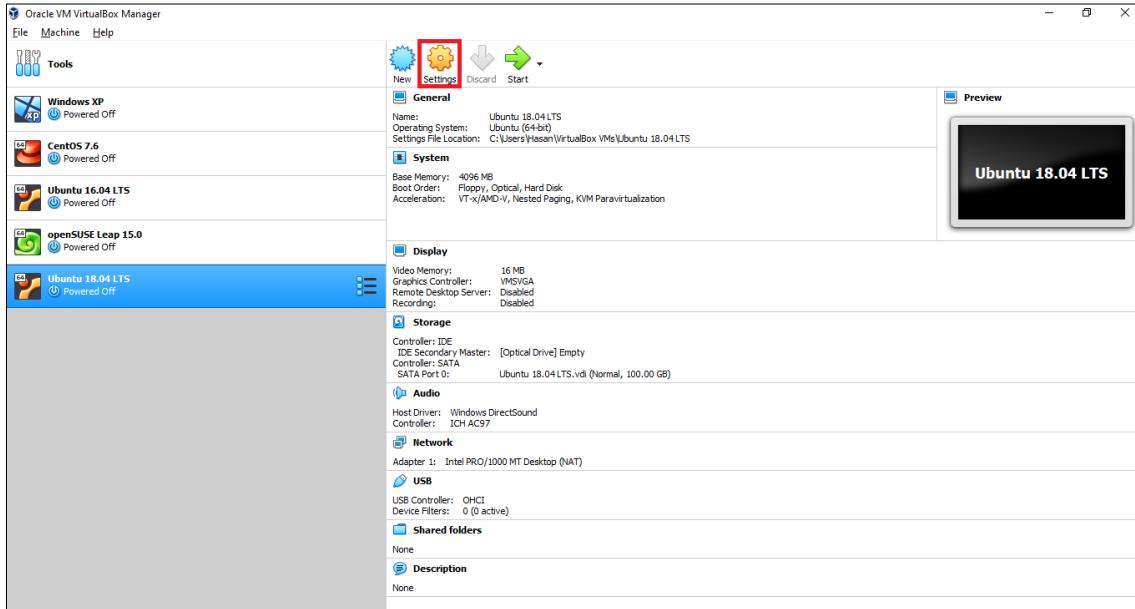


Figure 8: Select Settings

12. Select General → Advanced. Now, select the Shared Clipboard and Drag'n'Drop option to Bidirectional.

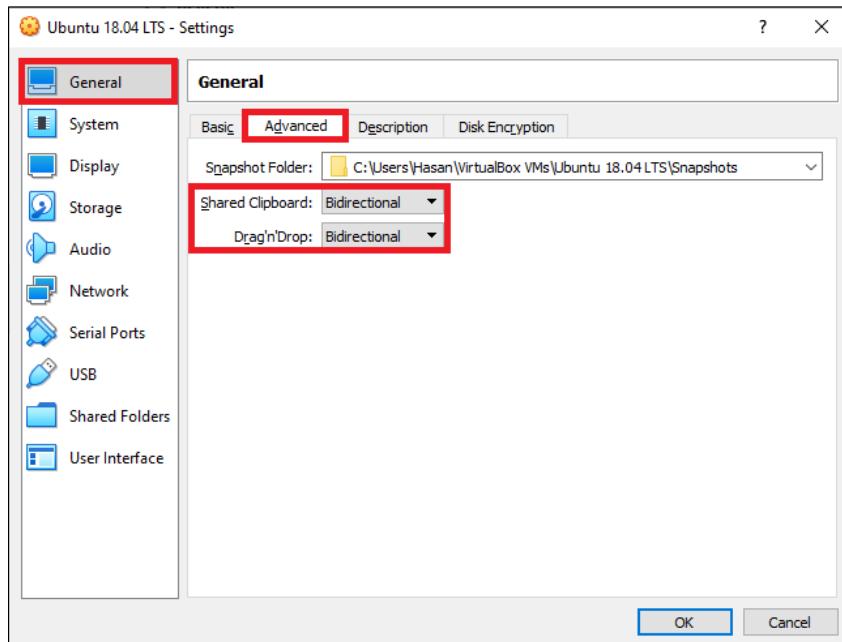


Figure 9: Advanced Settings

13. Select Settings to assign the image file of respective OS to VB. **Select Storage → Controller : IDE → Empty.** Now, in the Attributes tab, click on New Disk and provide the path of downloaded image file of Ubuntu OS. **Click OK.**

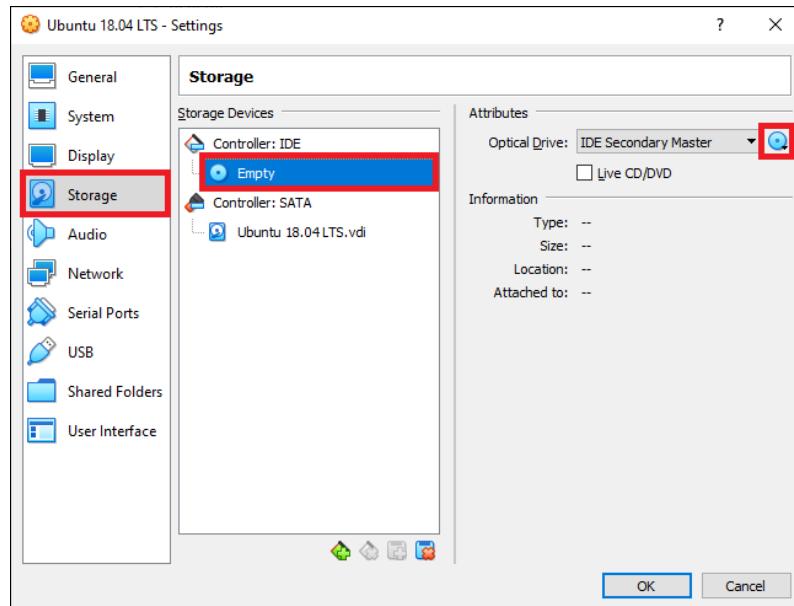


Figure 10: Providing Image file

14. Once everything is in place, it's time to boot that ISO and install Linux as a virtual operating system. **Click Start.**

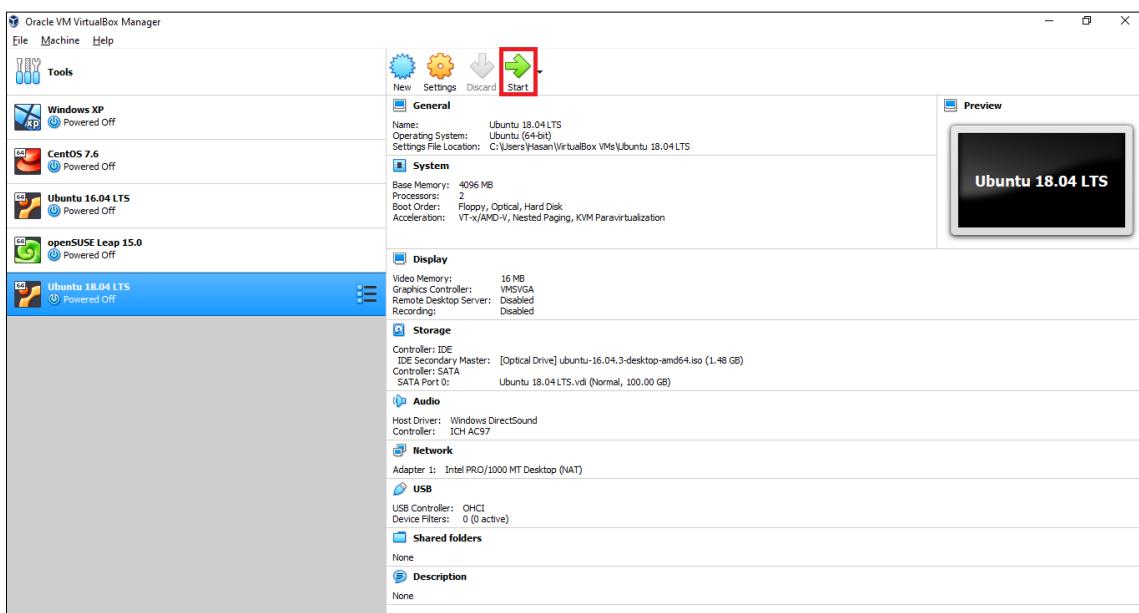


Figure 11: Starting Virtual OS

15. Virtual OS will boot into Linux Installation process. You should be presented with the option to install it. Click **Install Ubuntu**.

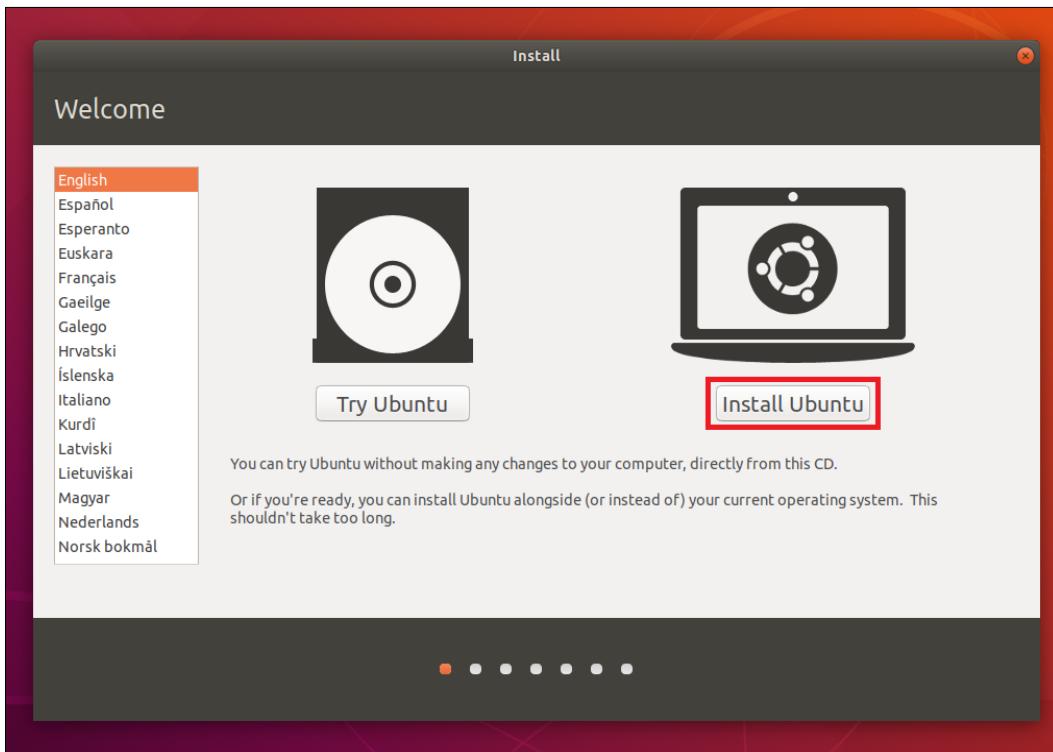


Figure 12: Installing Ubuntu

16. Continue with **Normal Installation**.

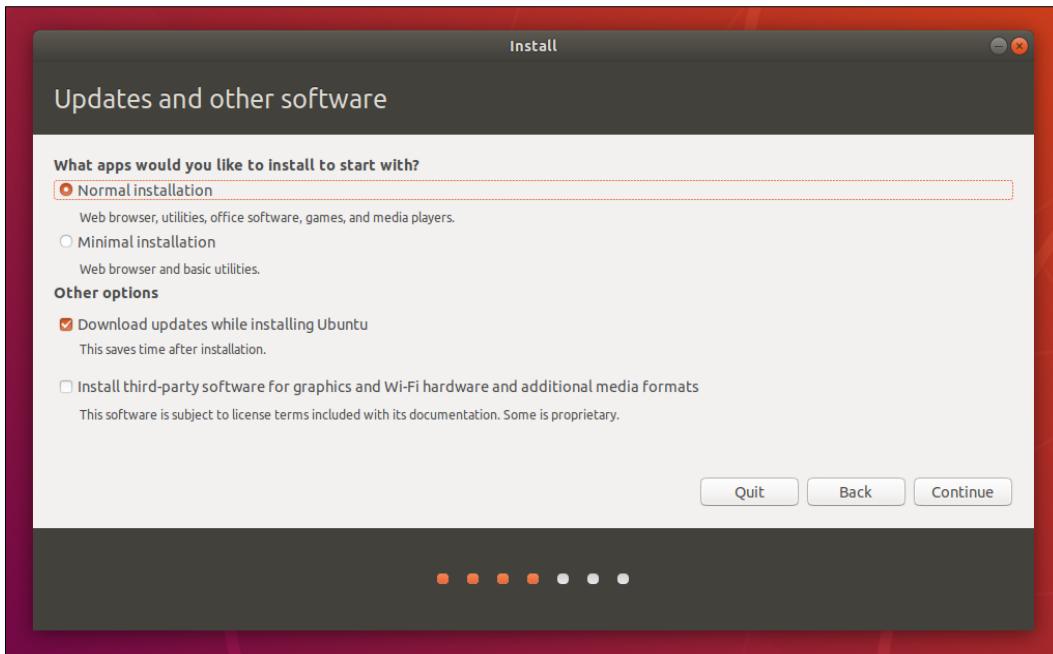


Figure 13: Update Screen

17. In Installation type screen, select **Erase disk** and install Ubuntu option.

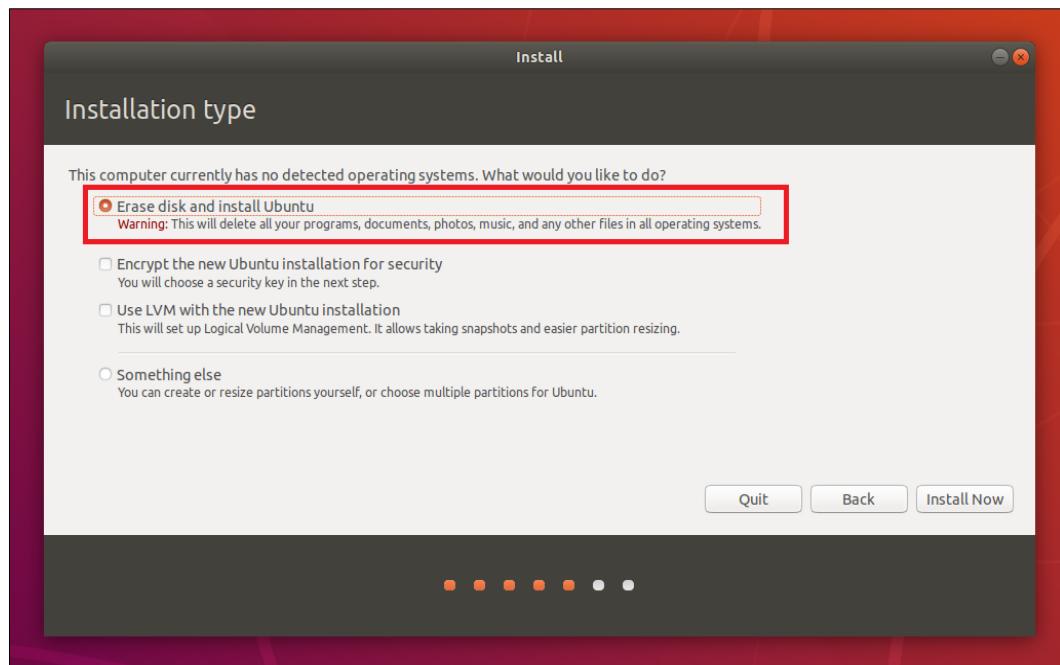


Figure 14: Installation type Screen

18. Select **Continue**.

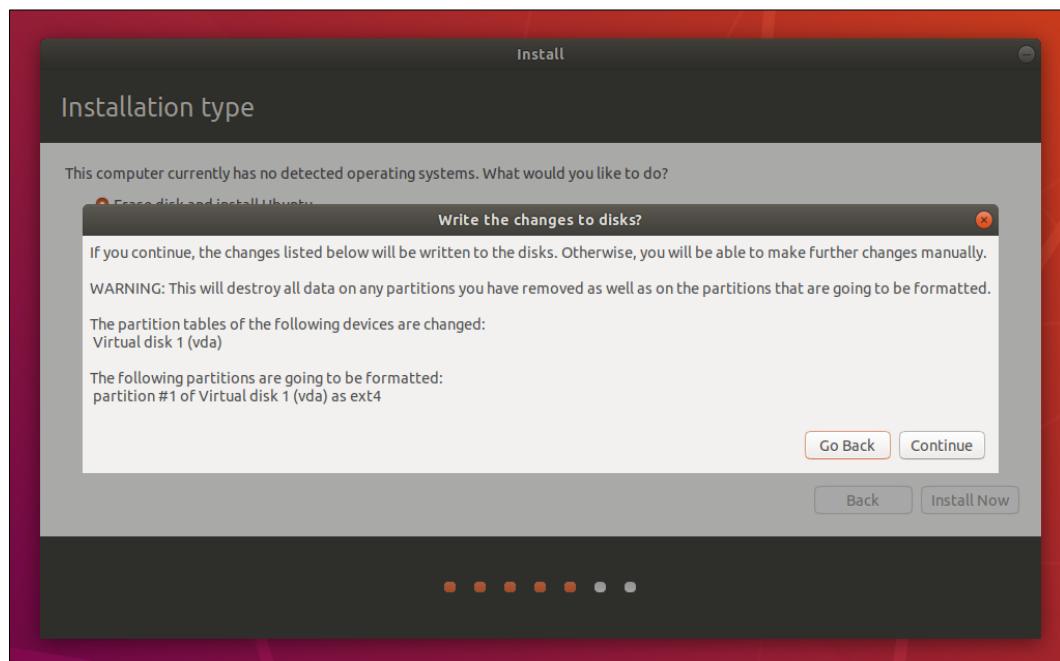


Figure 15: Confirmation Screen

19. Select your Current City and Continue. Select **Continue**.

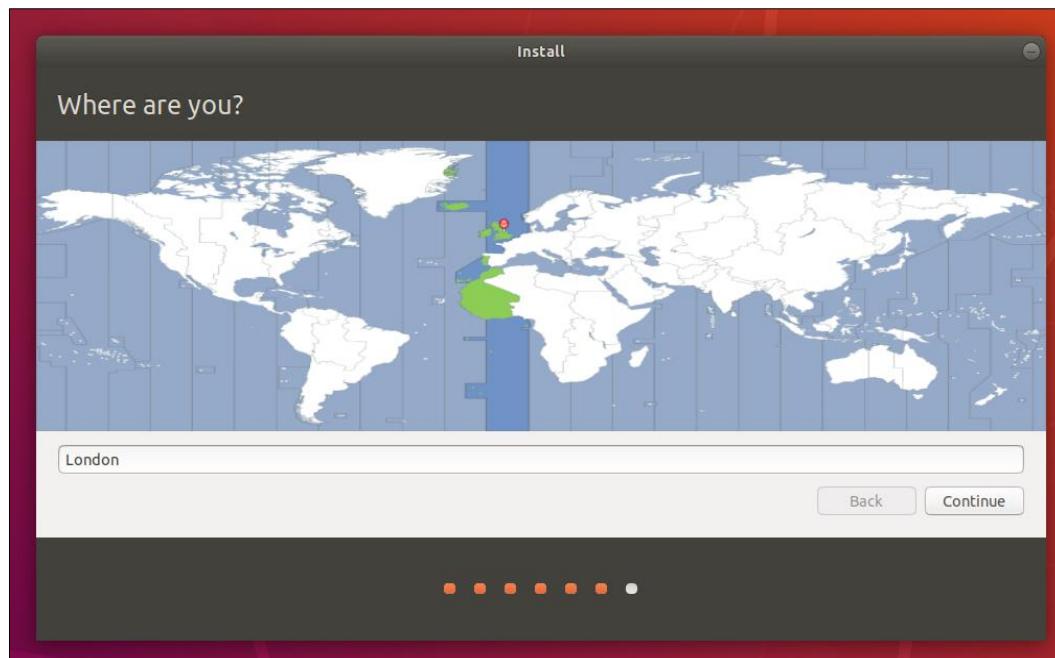


Figure 16: Location selection screen

20. Fill your Info and click **Continue**.

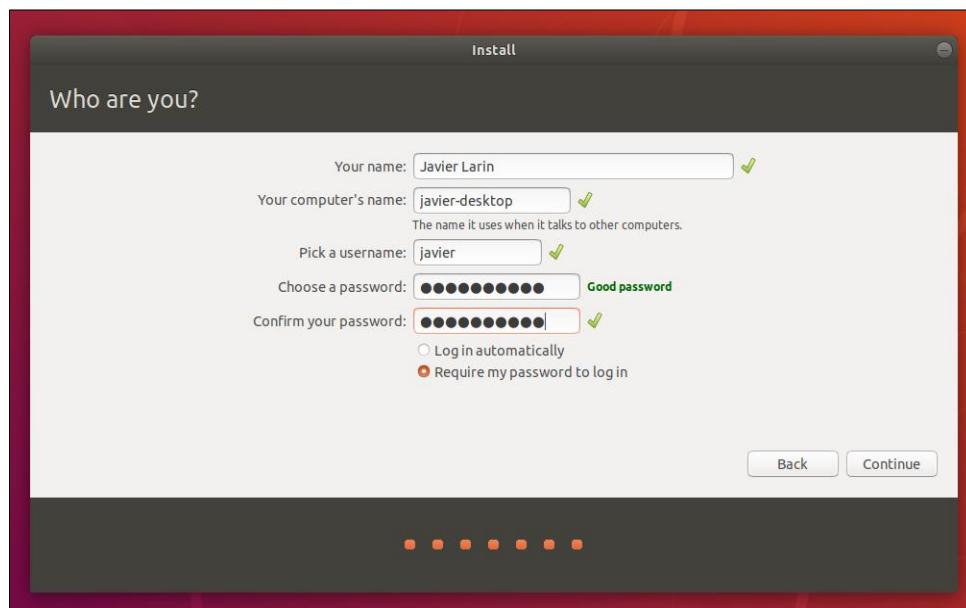


Figure 17: Intro Screen

21. Installation will Continue.

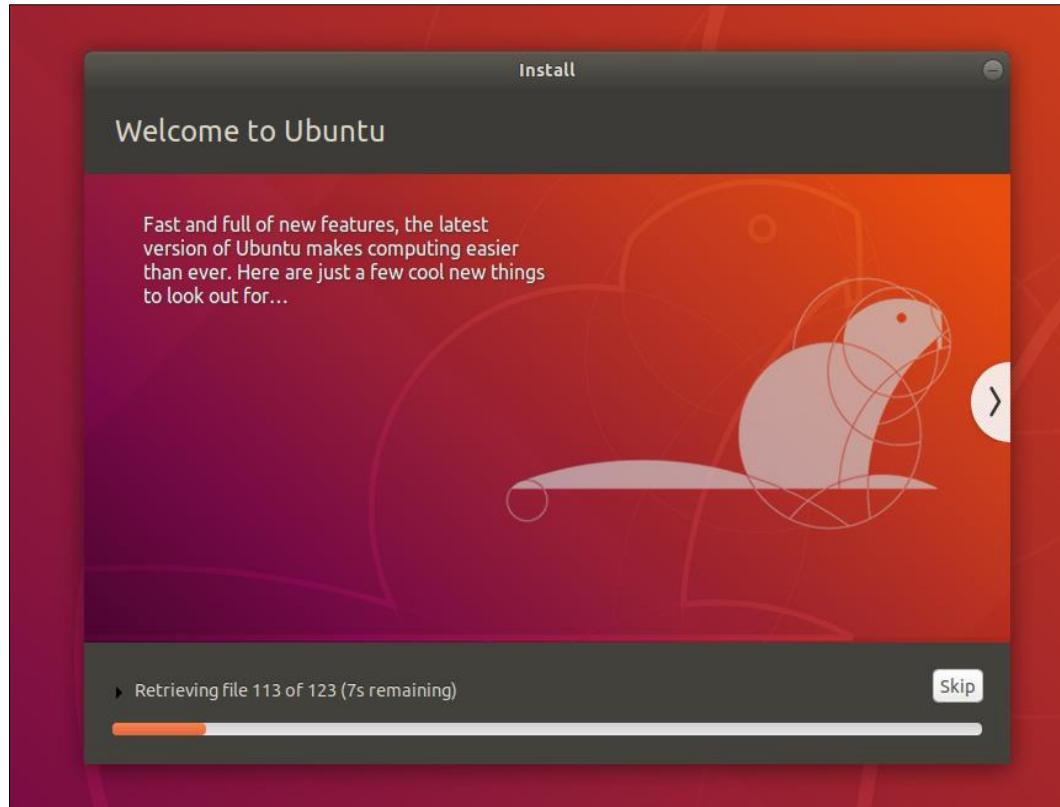


Figure 18: Installation Screen

22. Installation is Complete. Click **Restart Now**.

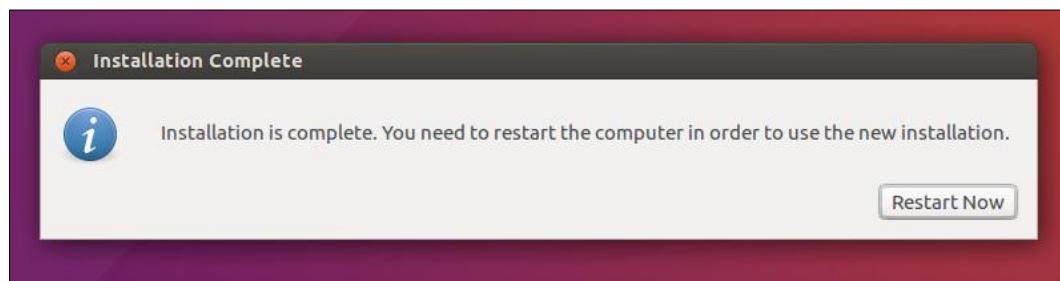


Figure 19: Complete Installation Screen

# aXonica-UI Application

---

## In this Chapter

We will learn about:

- aXonica graphical user interface
  - Benefits for users unaware of UNIX language
- 

## Introduction

The users who have a difficulty in understanding shell script don't have to learn anything about it because aXonica-UI provides an interactive interface for the users. It is basically developed using Qt programming that runs shell file in the background.

## Configuration of aXonica-UI Application

aXonica-UI is downloaded along with the shell file of aXonica. After installation, a desktop icon is created like this:

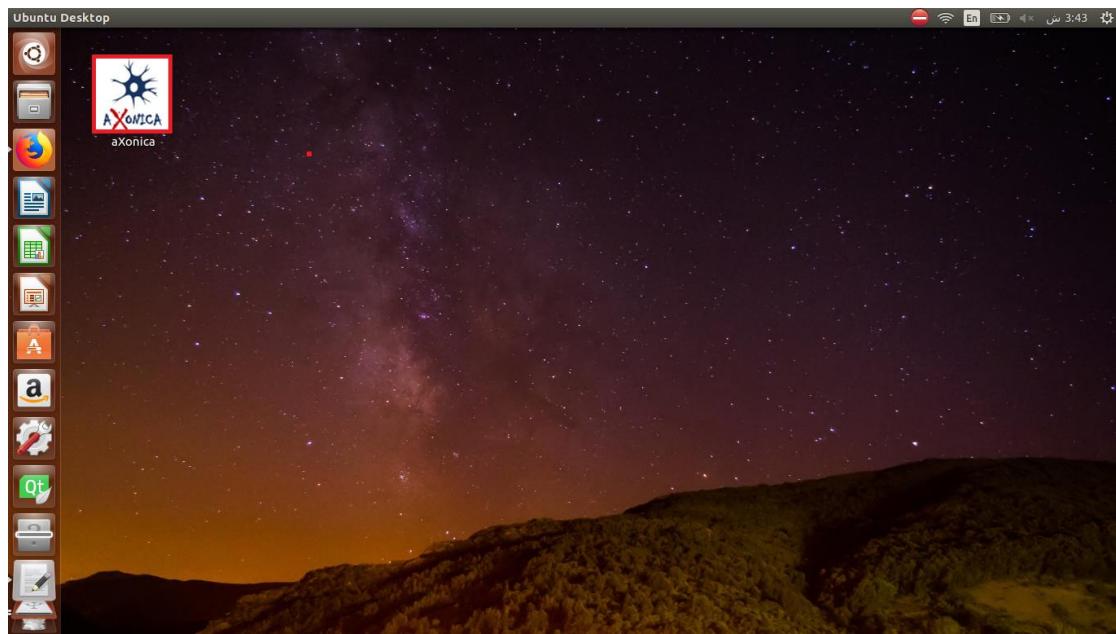


Figure 1: aXonica-UI desktop icon

1. Double click the icon. The following screen appears:

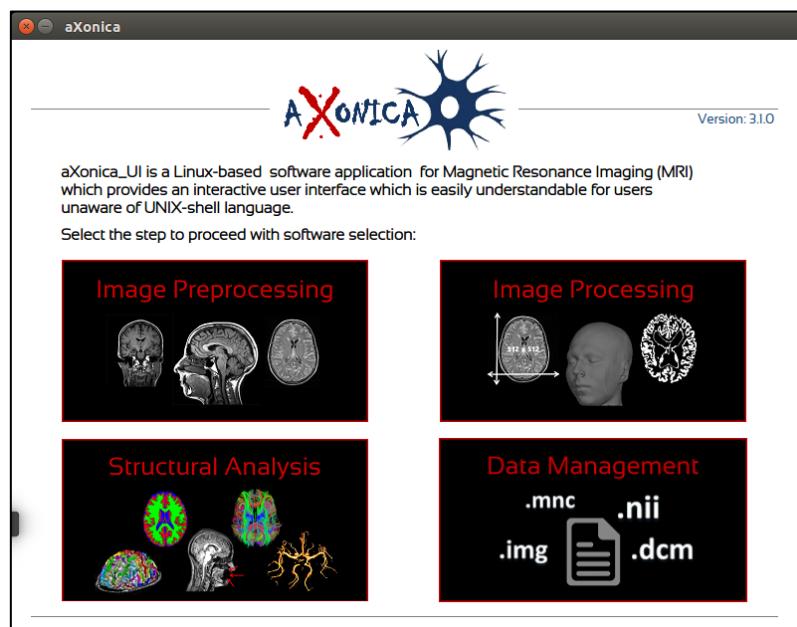


Figure 2: Main portfolio screen

2. Click on the category you wish to select (say Image Preprocessing). The following screen appears:

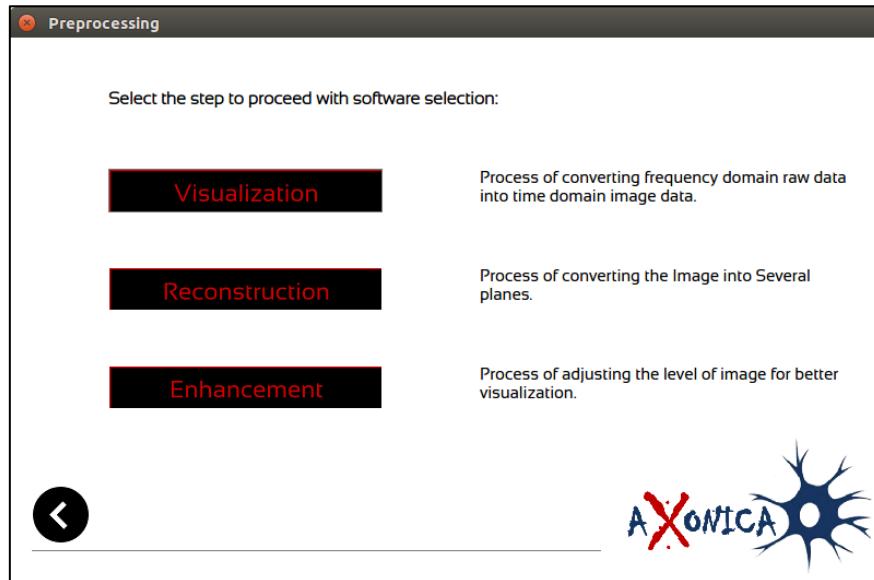


Figure 3: Sub portfolio screen

**NOTE** To go to the previous screen click 

3. After selecting the required category, the following screen appears:

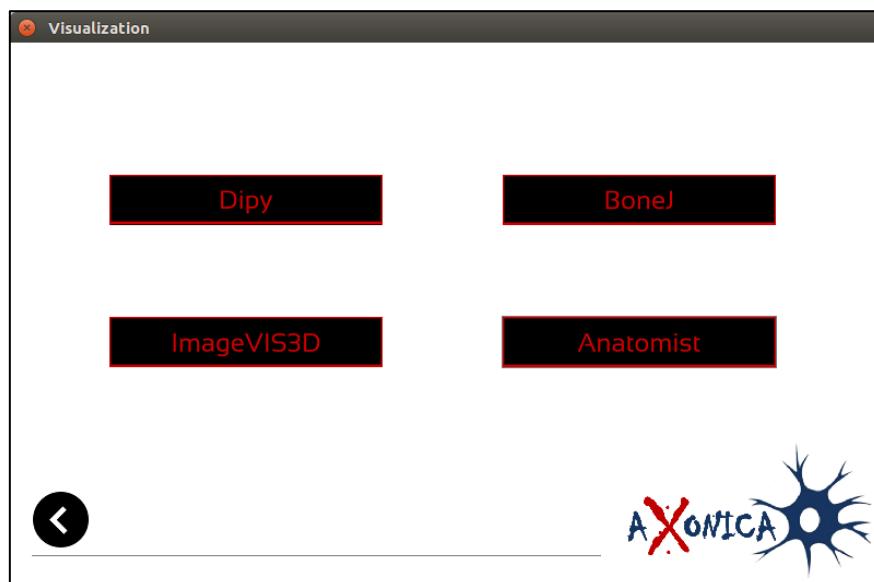


Figure 4: Tool selection screen

4. Click on your required software to run.

**NOTE** The buttons of the tools that were not downloaded by the user are disabled.

The procedure of running tools of other categories also proceeds in the same manner.

# Data Acquisition and Software Pipeline

---

## In this Chapter

We will learn about:

- Image Data Acquisition
  - MRI network diagram in hospitals
  - Introduction to software pipeline
- 

## Introduction

Medical Imaging Technology is identified globally by some major diagnostic imaging device manufacturers. Honorable mentions contain Fujifilm Holdings, Siemens Healthcare, Toshiba Medical Systems Corporation and Canon Medical Systems Corporation.

MRI machine is installed and distributed with some extra workstation for functionality i.e.

- Scan Console
- Display Console
- Reconstruction Box
- Quality Control (Calibration Phantoms)

## Image Data Acquisition from MRI

MRI installation network diagram is given as follows:

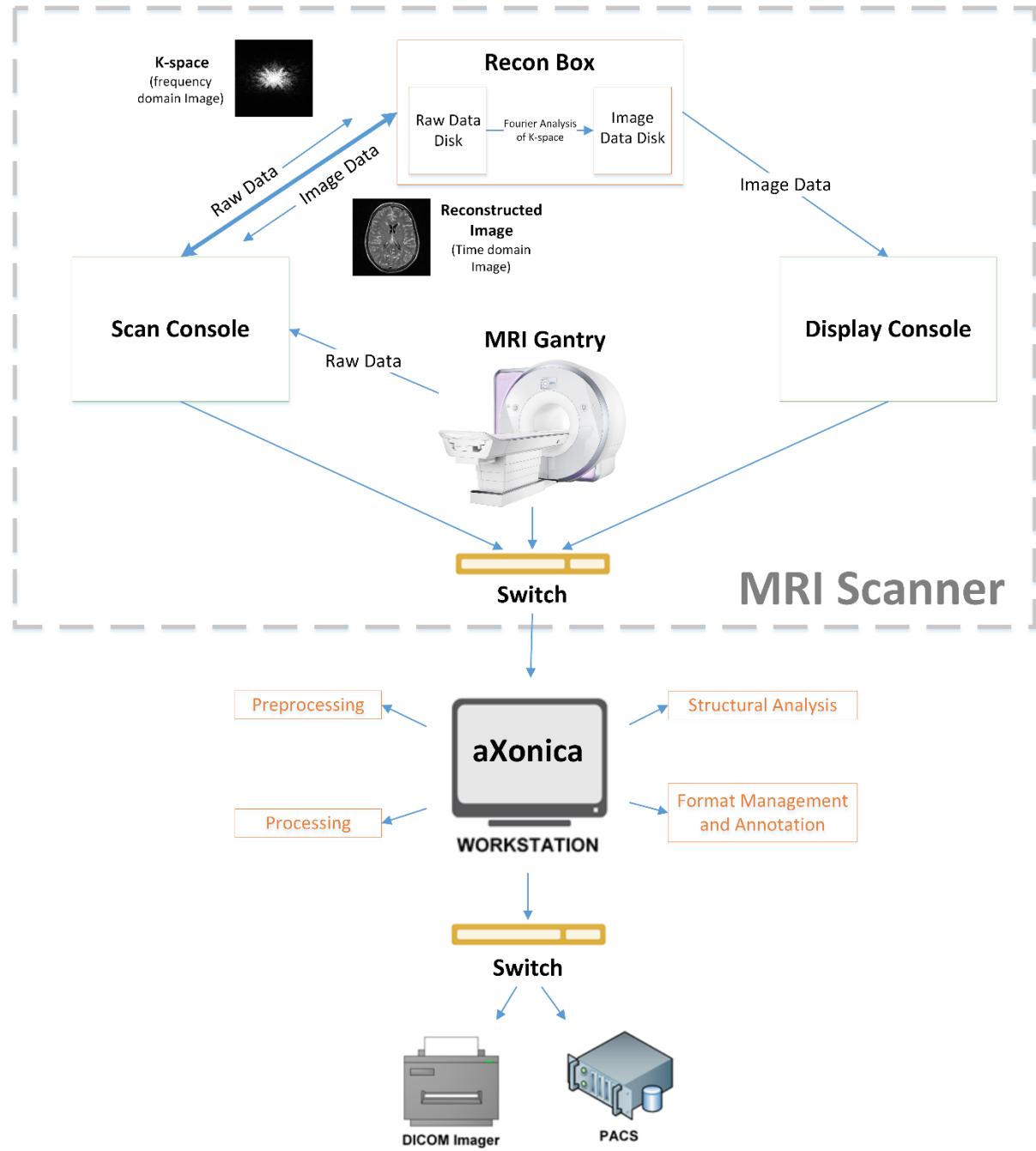


Figure 1: MRI network diagram

**STEP 1:**

In the first step, the Scan Console is given commands by the user about which part of the human body is about to be analyzed.

**STEP 2:**


These commands are transferred to the MRI and patient data is acquired. This raw data is in K-space (frequency domain).

**STEP 3:**

The data is then sent to the Scan Console which further sends this data to the reconstruction box. This raw data is stored in Raw Data Disk (RDD).

**STEP 4:**

Inside the reconstruction box, Fourier analysis is performed on the data and thus it is converted into time domain.

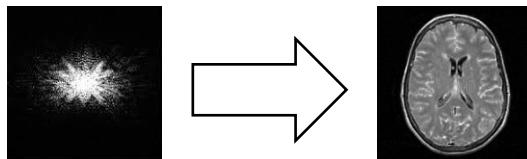

The output, known as Image Data received from the transform is saved in Image Data Disk (IDD).

**STEP 4:**

The image data acquired from IDD is then transferred to Display Console for further processing.

**STEP 6:**

Now the data present in Display Console is further transferred to aXonica workstation which has all freeware tools installed in it. Any operation can be performed according to user's request.

## Software Image Pipeline

MRI software base is divided into several steps. These steps are visualized using a pipeline so that the tools can be accessed according to user's needs. aXonica contains several tools from each category.

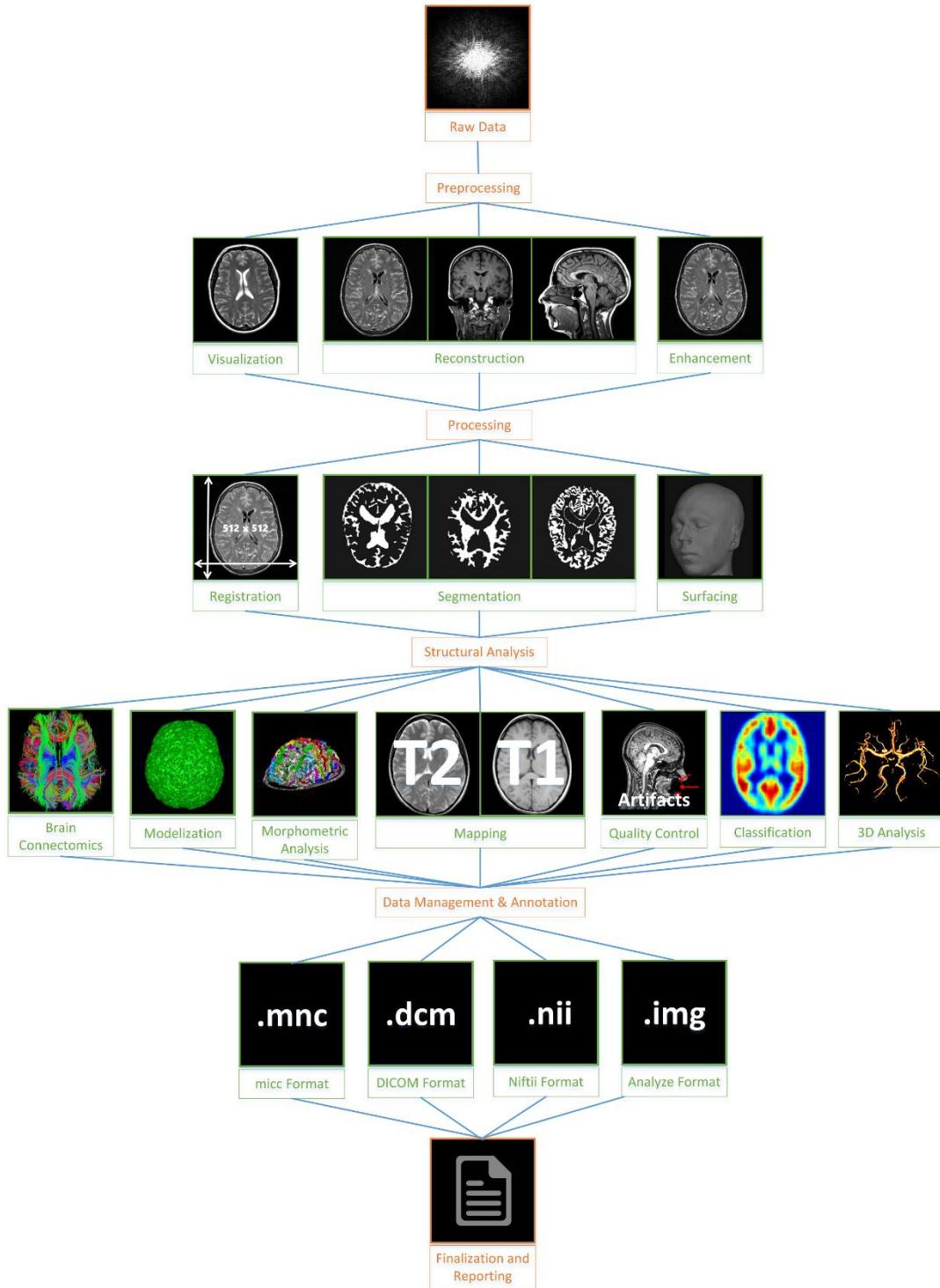


Figure 2: MRI Pipeline Visualized

# Tutorial Datasets

## Downloading aXonica Tutorial Datasets

aXonica tutorial datasets is available at its website:

<https://github.com/hasaniqbal777/aXonica-bin>

1. aXonica can also be downloaded from the following command through git:

```
$ git clone https://github.com/hasaniqbal777/aXonica-bin
```

2. Now run the following commands on terminal:

```
$ chmod +x aXonica_tutorial-datasets_enUS  
$ sudo bash aXonica_tutorial-datasets_enUS
```

3. Download wizard of tutorial datasets will start.

**NOTE** Download tutorial datasets *require a proper internet connection to proceed, otherwise the installation terminates.*

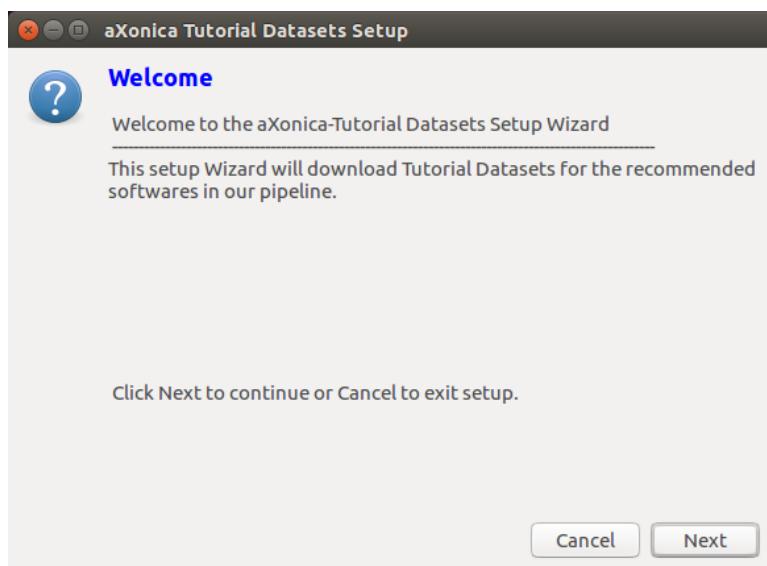


Figure 1: Welcome screen

4. Click **Next** to proceed and confirm the Internet connection.

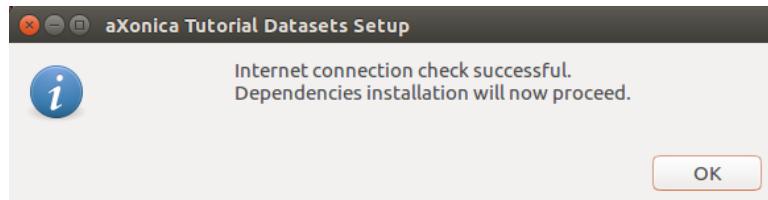


Figure 2: Internet check successful screen

5. Installation is **terminated** if there is **no internet**.

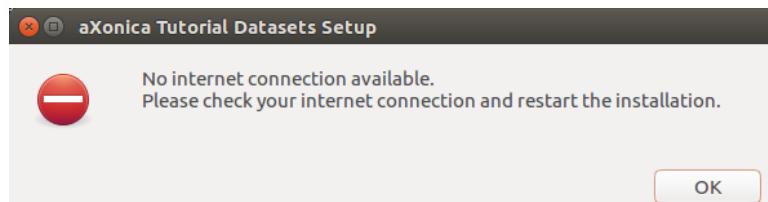


Figure 3: Internet check failed screen

6. Select the tutorial datasets you want to download.

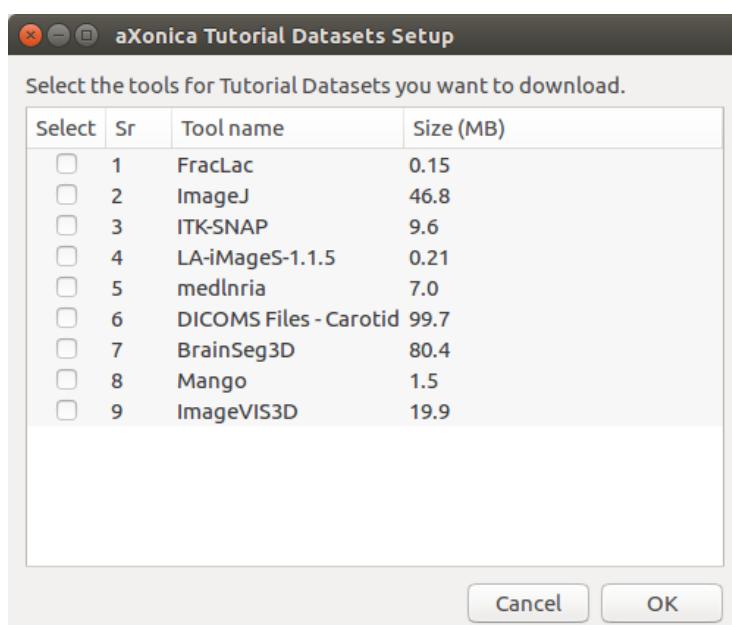


Figure 4: Tutorial Datasets selection screen

7. Click **Next** to proceed with the Download.

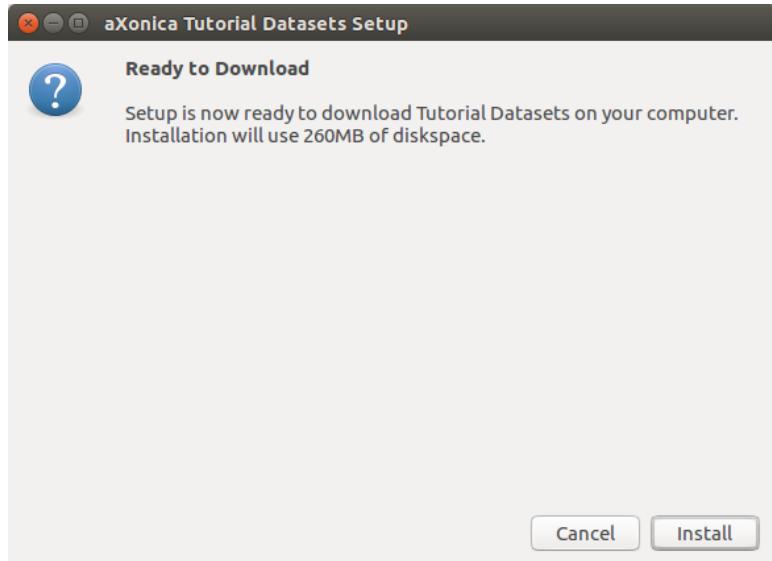


Figure 5: Ready to Download screen

8. Tools downloading will **continue**.


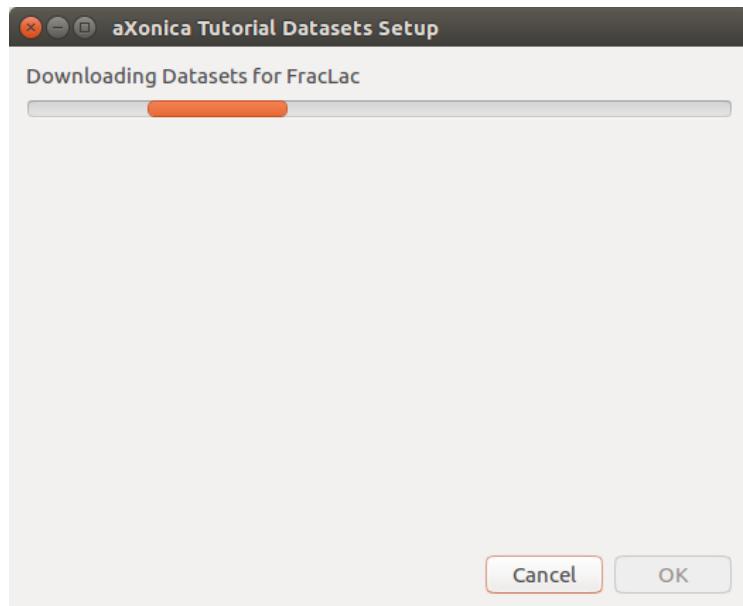


Figure 6: Download progress screen

9. Download of Tutorial Datasets is now finished. Click **Finish** to use the Tools.

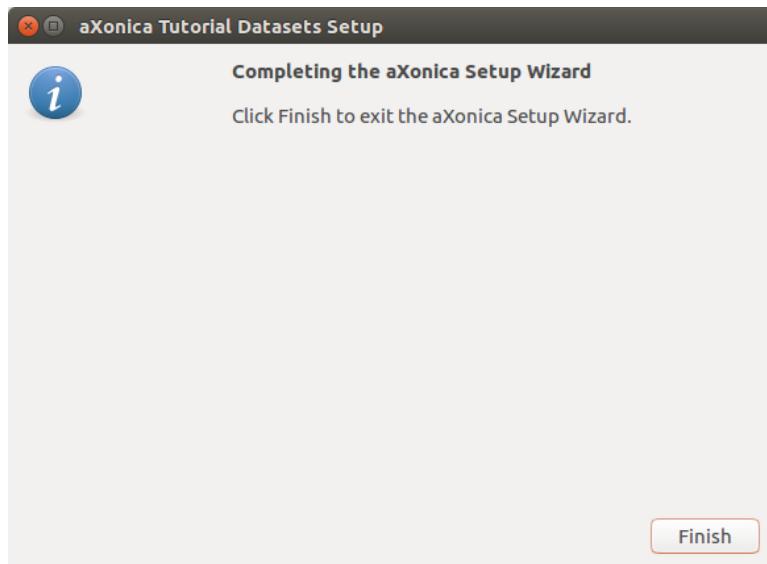

*Figure 7: Finalizing Installation screen*

# Image Preprocessing

---

## In this Chapter

We will learn about:

- Image Preprocessing
  - Steps for Preprocessing
  - Recommended tools for Preprocessing
- 

## Introduction

Pre-processing concerns the improvement of the Image Data and it enhances some image features that are important for further processing.

Image Preprocessing has three significant steps:

- Image Visualization
- Image Enhancement
- 3D Image Reconstruction

## Image Visualization

Raw data is acquired from MRI Scanner. The images acquired are in the frequency domain. Specific Fourier analysis converts the “Raw Data” into “Image Data” which is in time domain. It is further pre-processed and transferred for imaging.

**Recommended tool** *ImageVis3D*

### ImageVis3D [1]

**An architecture for Large Scale Volume Rendering:** Provides domain specific visualization capabilities.

1. Open terminal and run these commands:

```
$ cd /usr/local/ImageVis3D-3.1.0/
$ ./ImageVis3D
```

2. Select “Open data set from file”


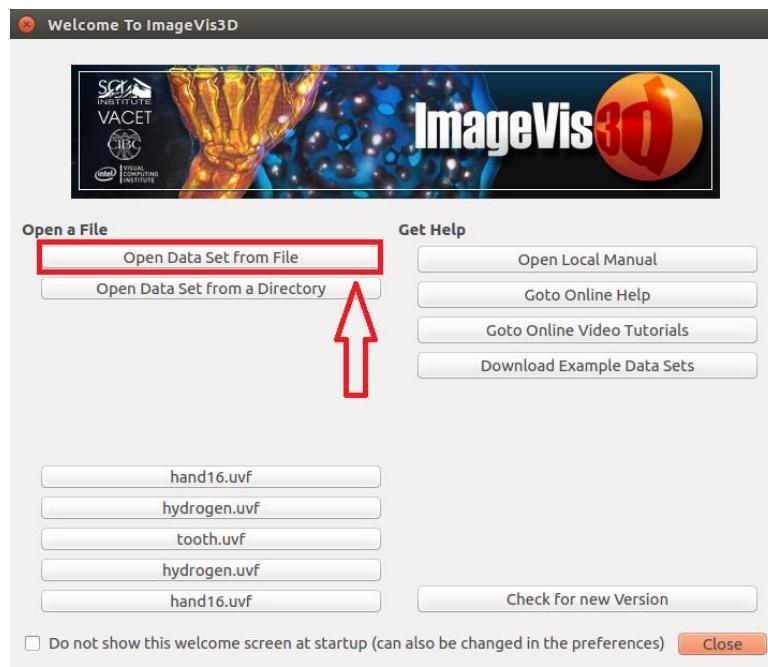


Figure 1: ImageVis3D Welcome screen

**3. Load the downloaded dataset**

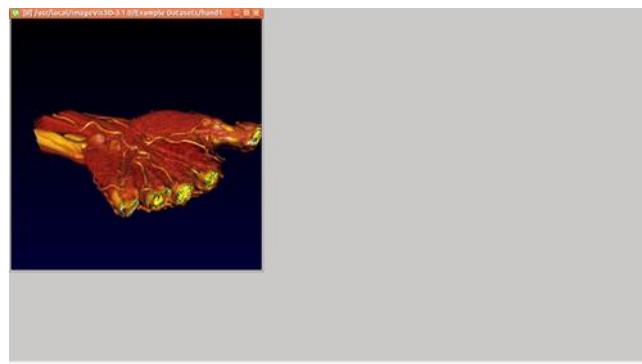

*Figure 2: Visualized Image*

**4. Further Analyze the image.**

## Image Enhancement

Brightness of specific parts of body (e.g. Brain, Spine etc.) is adjusted so that the results are more suitable for display and further image analysis.

**Recommended tool** **AMIDE**

### AMIDE [2]

**A Free Software Tool for Multimodality Medical Image Analysis:** Displays and analyzes multimodality volumetric medical images.

1. Run the following commands in the terminal:

```
$ export UBUNTU_MENUPROXY=0
$ amide
```

2. AMIDE "New Study" screen opens:

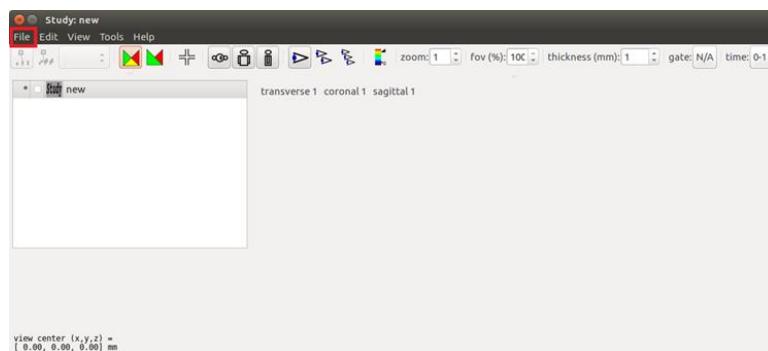


Figure 3: AMIDE Welcome screen

3. File → Open Study

4. Load the downloaded tutorial dataset file: *m2862-small.xif*

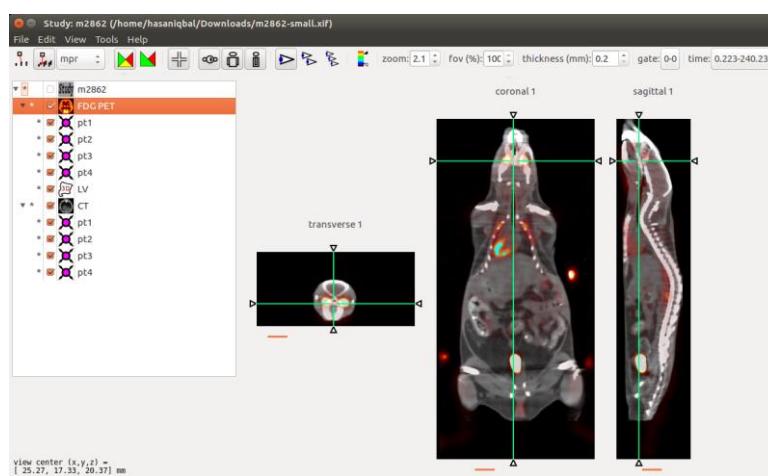


Figure 4: Enhanced Images

5. Further Analyze the image.

## 3D Image Reconstruction

Visualized images are processed and several images are aligned and converted to separate Axial, Coronal and Sagittal planes of the visualized part of body.

**Recommended tool** *ITK-SNAP*

### ITK-SNAP [3]

**Allows users to segment structures in 3D medical images.** ITK-SNAP provides semi-automatic segmentation using active contour methods, as well as manual delineation and image navigation.

1. Run the following commands in the terminal:

```
$ cd /usr/local/itksnap-3.6.0-20170401-Linux-x86_64
$ cd bin/
$ ./itksnap
```

2. File → Open main image → Browse

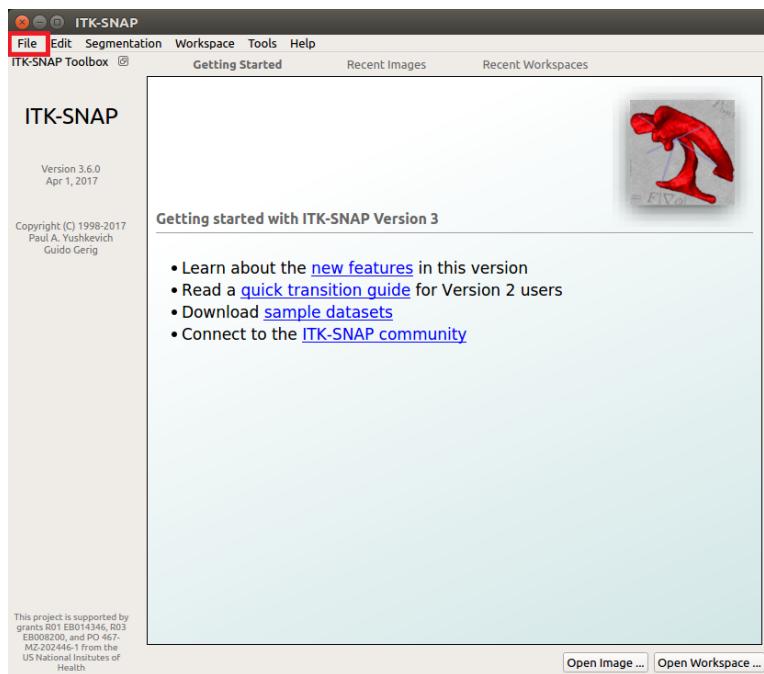


Figure 5: ITK-SNAP Welcome screen

### 3. Load data set from file:

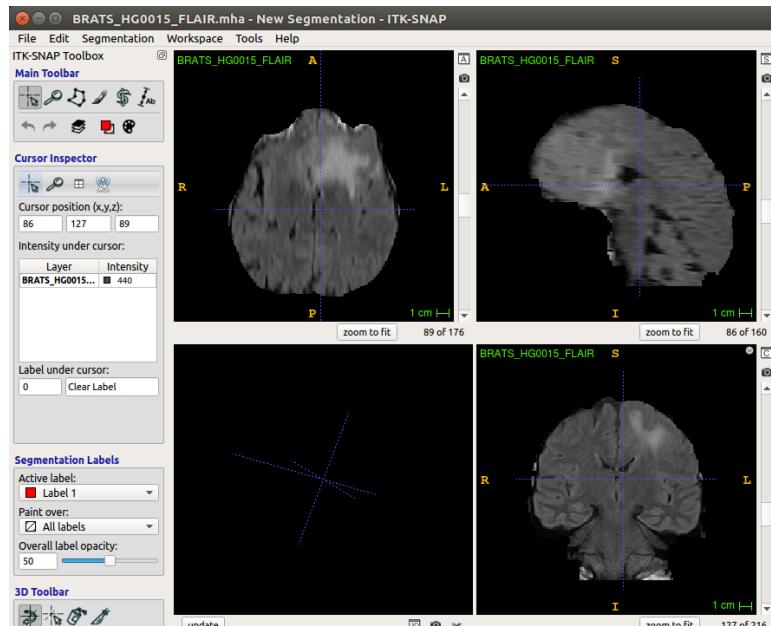


Figure 6: 3D Reconstructed Image

### 4. Further Analyze the image

# Image Processing

---

## In this Chapter

We will learn about:

- Image Processing
  - Steps for Processing
  - Recommended tools for Processing
- 

## Introduction

Image processing refers to processing or altering an existing image in a desired manner. Image Processing has three significant steps:

- Image Segmentation
- Image Registration
- Image Surfacing

## Image Segmentation

It is a process of partitioning a digital image into multiple segments which are more meaningful and easier to analyze. For Example: A Brain MRI image is segmented into several images for visualization and analyses of gray matter, white matter and cerebrospinal fluids.

**Recommended tool** *BrainSeg3D*

### BrainSeg3D [4]

Provides a free volume (3D image) viewer and segmentation tool. BrainSeg3D is a graphic application that make segmentation of volumes more accurate by providing tools for semi-automated segmentation combined with a user friendly graphic interface.

1. Run these commands in the terminal:

```
$ cd /usr/local/brainseg3d
$ ./BrainSeg3D
```

2. Select 'Quick open file'

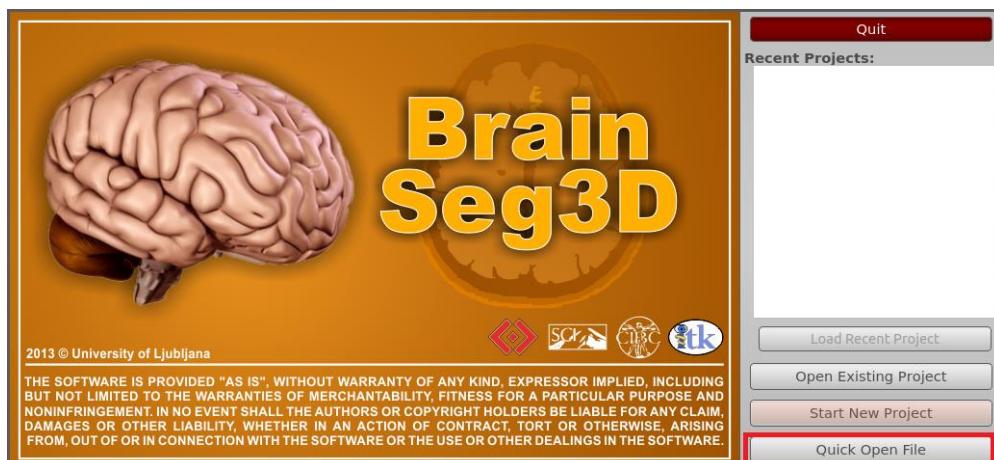


Figure 1: BrainSeg3D Welcome Screen

3. Load Tutorial data sets

#### 4. Select all Images and Import

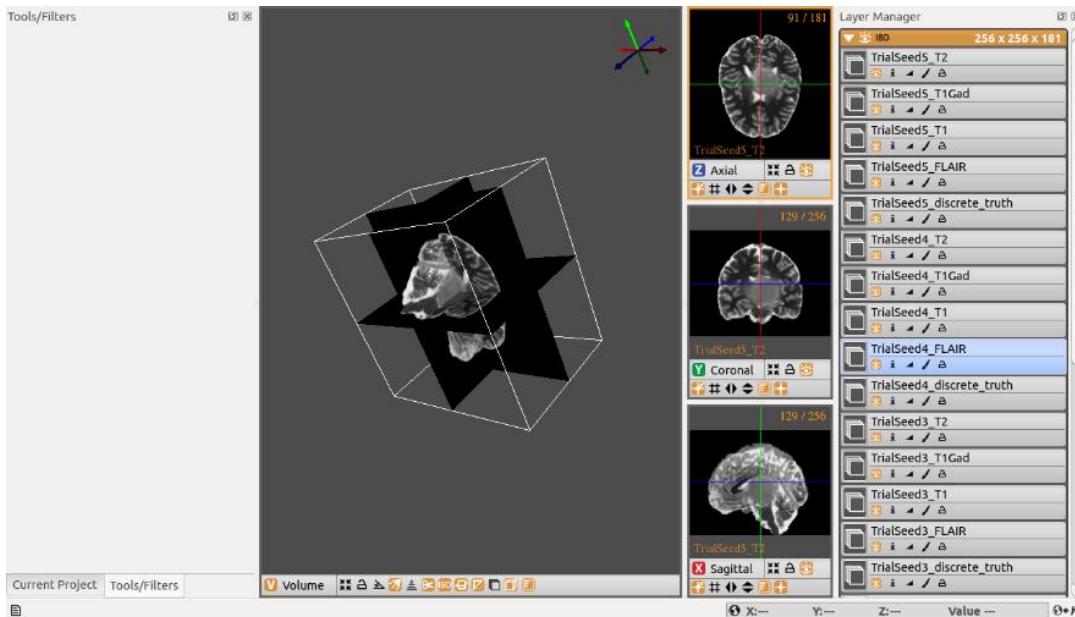


Figure 2: Segmented Image

#### 5. Further Analyze the image

## Image Registration

Process of transforming different sets of data into one coordinate system.

Standard Co-ordinate system for DICOM images is 512 x 512-pixel resolution. This provides a standard base for patient imaging. It is very helpful in reporting of the patient to have a standard pixel ratio.

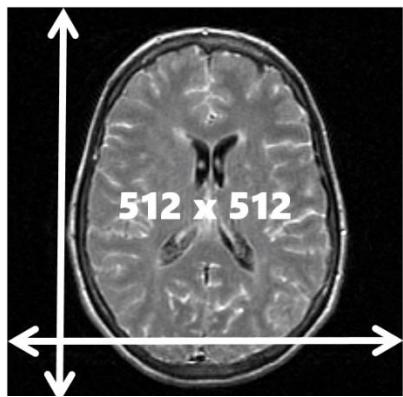


Figure 3: Image Registration

**Recommended tool** *Mango*

## Mango [5]

A complete data analysis pipeline that provides statistical confidence estimates for interactions and corrects for major sources of bias including differential peak enrichment and genomic proximity.

1. Run the following commands in the terminal:

```
> cd /usr/local/Mango
> ./mango
```

2. Select **Open**.

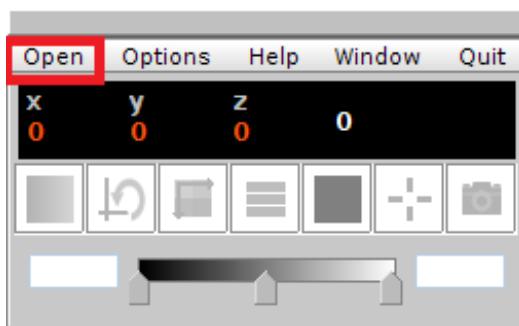


Figure 4: Mango Welcome screen

3. Open the Tutorial Dataset downloaded.

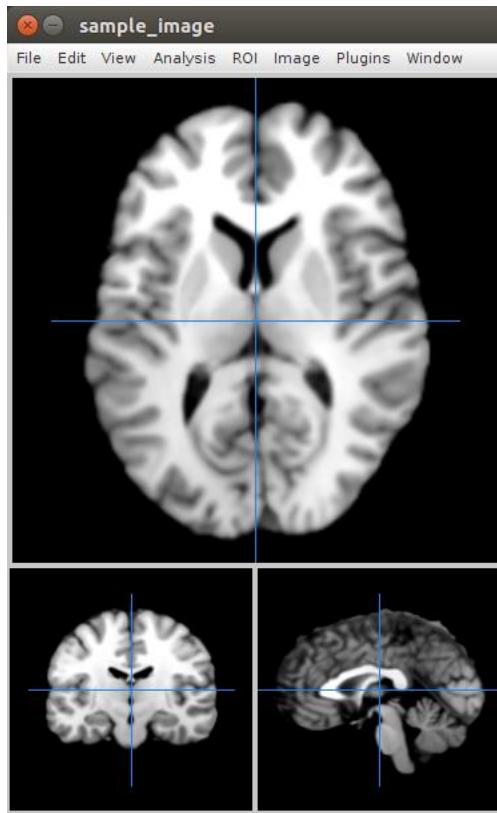


Figure 5: Registered Image

4. Further Analyze the image according to needs.

## Image Surfacing

Process in which a 3D view of the Image Data is created and visualized.

**Recommended tool** **BrainVISA**

### BrainVISA [6]

Neuroimaging Research Software Hosts heterogeneous tools dedicated to neuroimaging research.

1. Run the following commands in the terminal:

```
> cd /brainvisa-4.5.0
> ./anatomist
```

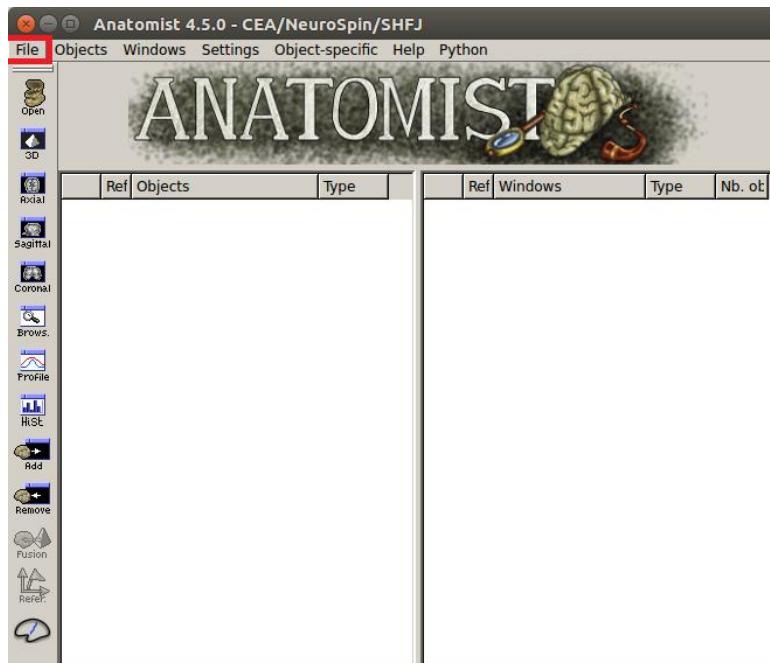


Figure 6: Anatomist Welcome screen

2. File → Open
3. Load the downloaded tutorial dataset file:
  - *T1 MRI: data\_for\_anatomist/subject01/subject01.nii*
  - *Activation map: data\_for\_anatomist/subject01/Video-T\_map.nii*
4. Select the 2 objects in the object list using Ctrl + left button
5. Click on the fusion button

6. A new window is **displayed** which allows to select some fusion parameters.

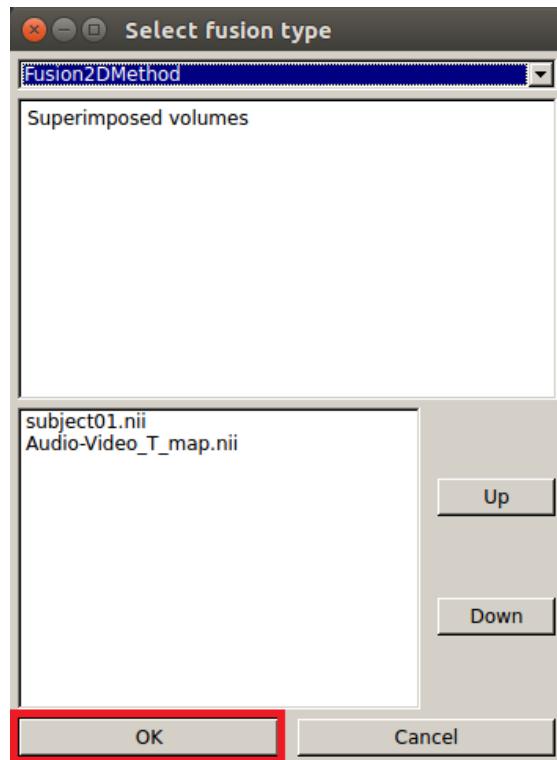


Figure 7: Fusion type selection screen

7. Click OK to create the fusion object

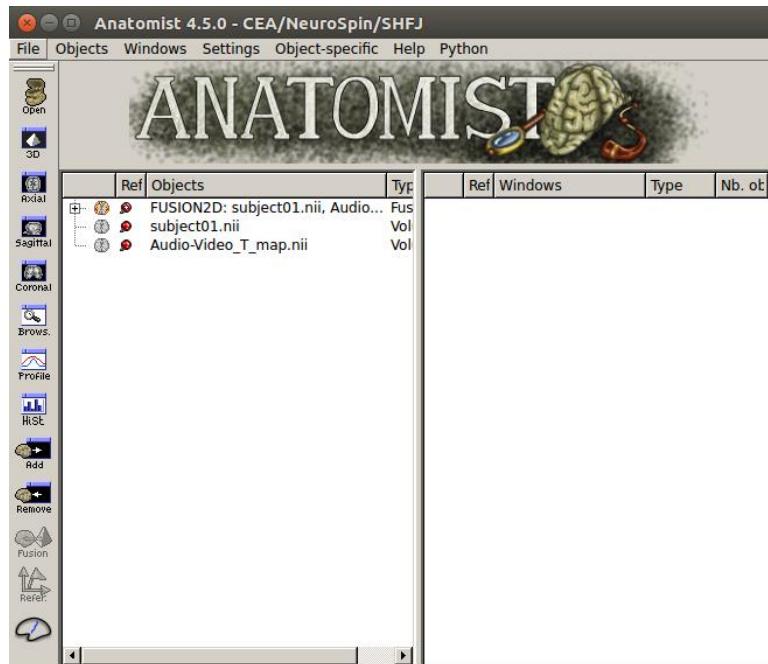


Figure 8: Fusion object screen

8. Open an axial window by clicking on 
9. Put the object (the fusion volume) into this window: **drag and drop** this object into the window. Following non registered image will appear

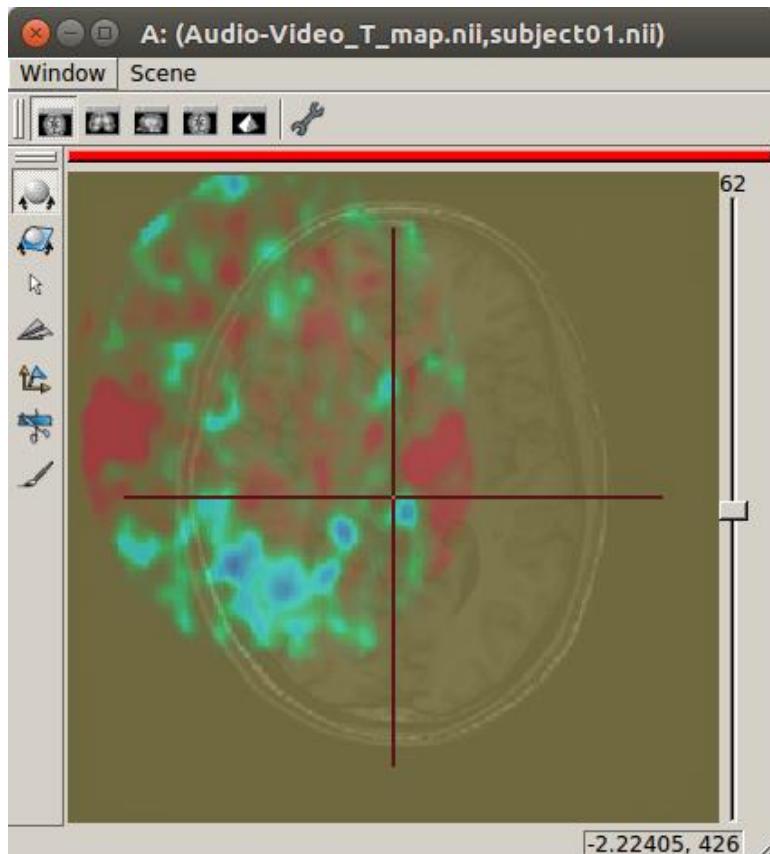


Figure 9: Fused mismatched image

10. For each volume, **right-click**.
11. Referential → Load information from file header

**NOTE** *The option Referential → Load information from file header extracts information about referential and transformations which are stored in the image files i.e. automatic registration.*

12. Look at the fusion, the 2 images are now well superimposed.

**NOTE** *Manual registration can also be done but in fact the human eye cannot drive a registration as well as a specific algorithm. For instance, images may seem aligned in an axial slice, but contain some drifts in sagittal and coronal orientations.*

# Structural Analysis

---

## In this Chapter

We will learn about:

- Structural Analysis
  - Steps for Structural Analysis
  - Recommended Tools for Structural Analysis of Images
- 

## Introduction


Structural analysis of an image refers to the extraction of meaningful information from the image by means of several techniques.

Structural Analysis has seven significant steps:

- Brain Connectomics
- Modelization
- Morphometric Analysis
- Quality Control
- Classification
- 3D Image Analysis
- Image Mapping


## Brain Connectomics

Brain Connectomics is the production and study of connectomes. Connectomes are comprehensive maps of connections within an organism's nervous system. Brain Connectomics Analysis is divided into two types:

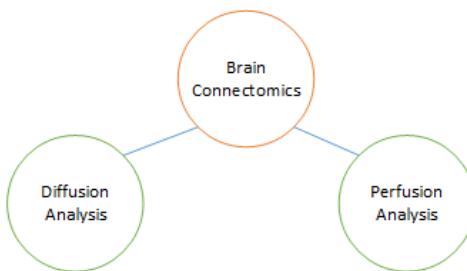

**Diffusion Analysis** is done by analyzing the diffusion of water molecules in tissues. Diffusion coefficient reduces for the cancerous tissues which is visualized.

**Perfusion Analysis** is done by analyzing the amount of blood taken up in certain areas of your brain. This can provide information on how your brain is functioning. Damaged areas have less oxygen and glucose use hence less blood supply which is visualized.

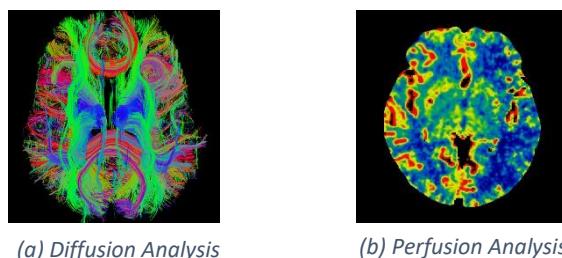

*Figure 10: Comparison Between Diffusion and Perfusion Analysis*

**Recommended tool** *medInria*

## medInria [7]

**Research and visualize medical images:** medInria is a platform for the diffusion of research software in medical imaging created by medInria teams.

1. Run the following commands on the terminal:

```
$ cd /usr/local/medinria-2.2.3-Linux-x86_64  
$ cd bin/  
$ ./medInria_launcher.sh
```

2. Select Diffusion

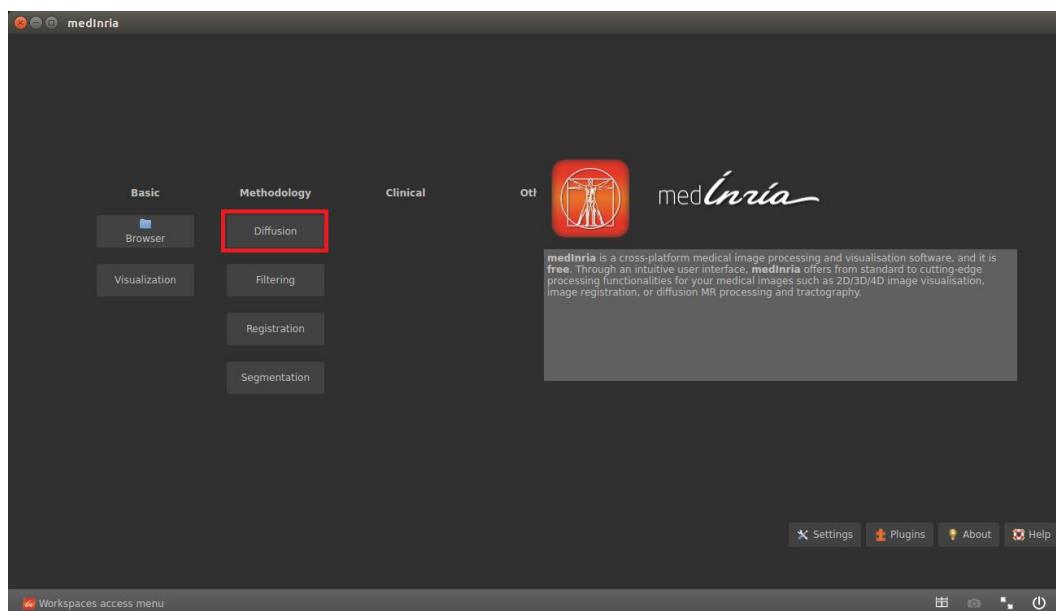


Figure 11: medInria Welcome Screen

5. Select "Open a file from your system".

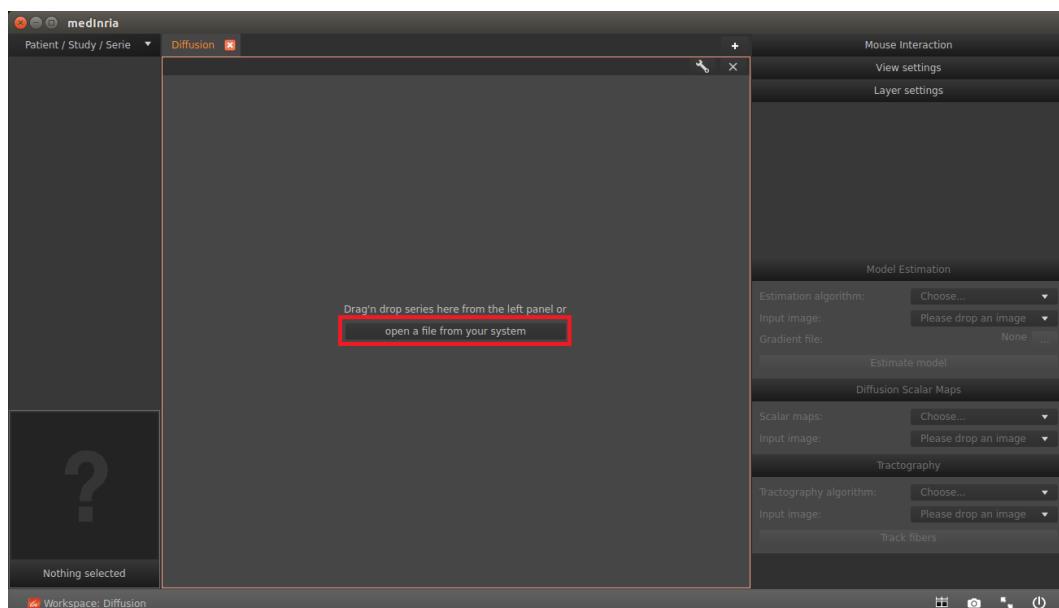


Figure 12: Load Study screen

5. Select the downloaded tutorial dataset file: ***xDTICUBE-4D.nii.gz***
6. Select the Estimation Algorithm as DTI estimation, in the Model Estimation tab.
7. Select the gradient file: ***gradient7.txt***



Figure 13: Model Estimation

8. Select Estimate model and wait for the diffusion analysis of the image

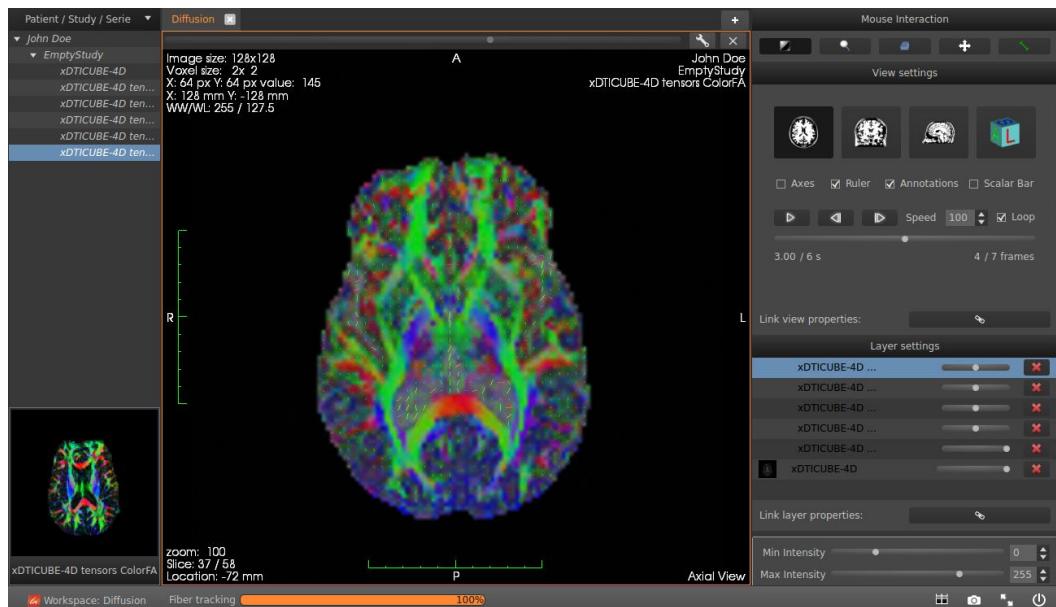


Figure 14: Diffusion Estimated Image

9. Further Analyze the Image according to needs.

## 3D Image Analysis

3D image analysis is the extraction of specific part of body from 3D visualized images which can be analyzed. For Example: 3D Analysis of Circle of Willis inside the brain.

Recommended tool **ImageJ**

### ImageJ [8]

**Biological Image Analysis:** Provides easy installation on arbitrary platforms and a simple user interface.

1. Run the following commands in the terminal:

```
$ cd /usr/local/ImageJ  
$ ./ImageJ
```

2. File → Open

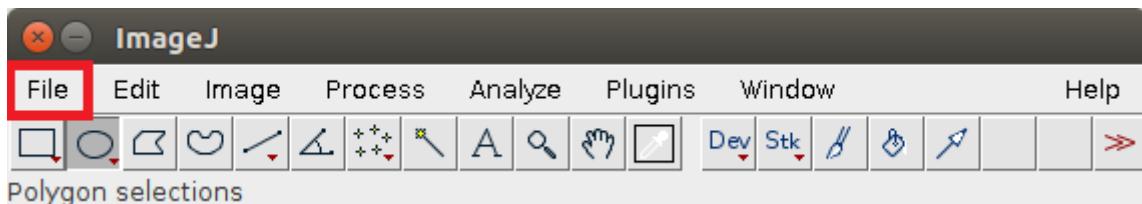


Figure 15: ImageJ Welcome Screen

3. Select the downloaded tutorial datasets: *t1-rendering.zip*

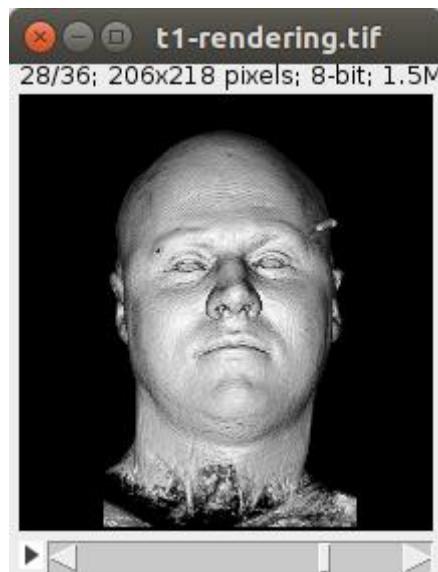


Figure 16: 3D Analyzed Image

4. Further Analyze the image according to needs

## Morphometric Analysis

It is visualization of the outer surface of the part of body under analysis.

**Recommended tool** *FracLac*

### FracLac [9]

**NOTE** *It is an ImageJ plugin, and it is executed from ImageJ.*

1. Run these commands in the terminal:

```
$ cd ~/Package/Plugins
```

**NOTE** *This is the resource folder of aXonica-master where the tools are downloaded.*

```
$ sudo cp Frac_Lac.jar /usr/local/ImageJ/plugins
```

**NOTE** *Provide the password for administrator.*

```
$ cd /usr/local/ImageJ
```

```
$ ./ImageJ
```

2. File → Open

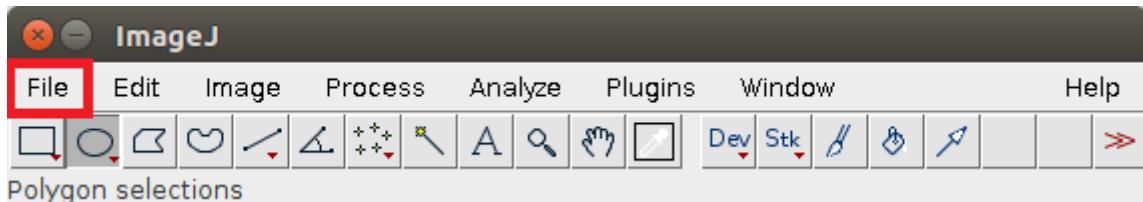


Figure 17: ImageJ Welcome Screen

3. Select the downloaded tutorial datasets: testimage.gif



Figure 18: Sample Image

### 3. Select "BC"

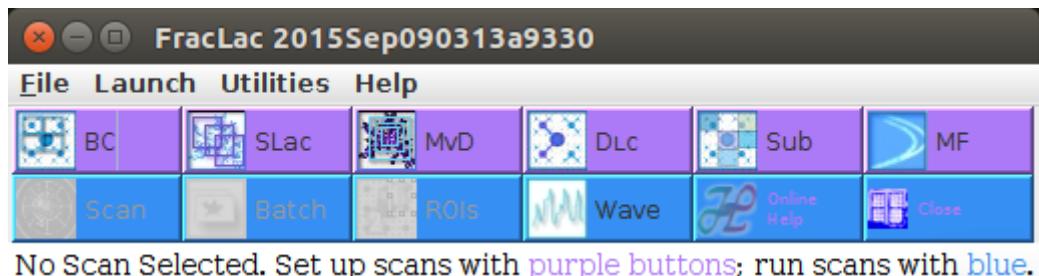


Figure 19: FracLac Welcome screen

### 4. Select OK



Figure 20: Mode Selection

### 5. Select OK

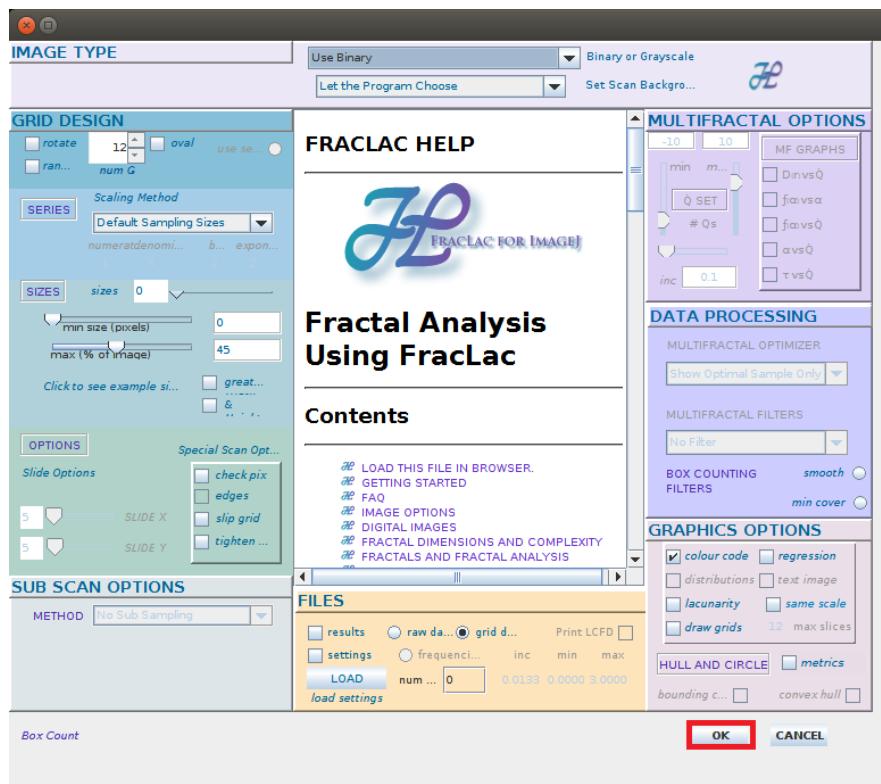


Figure 21: FracLac settings screen

## 6. Select OK

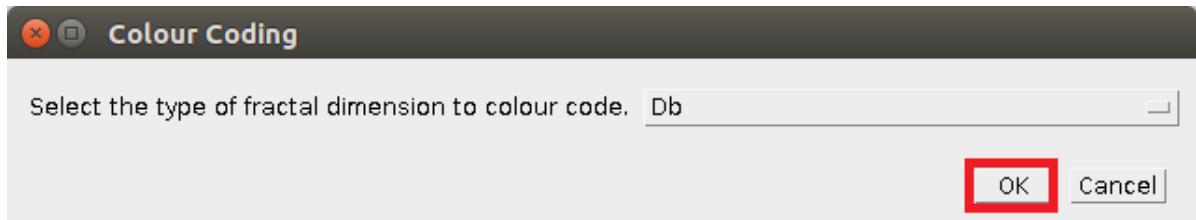


Figure 22 Colour Coding Selection screen

## 7. Select Scan

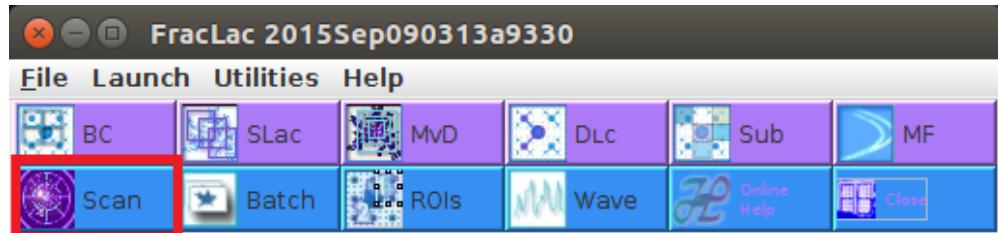


Figure: 23 Ready to scan

**NOTE** Wait for the image to be processed.

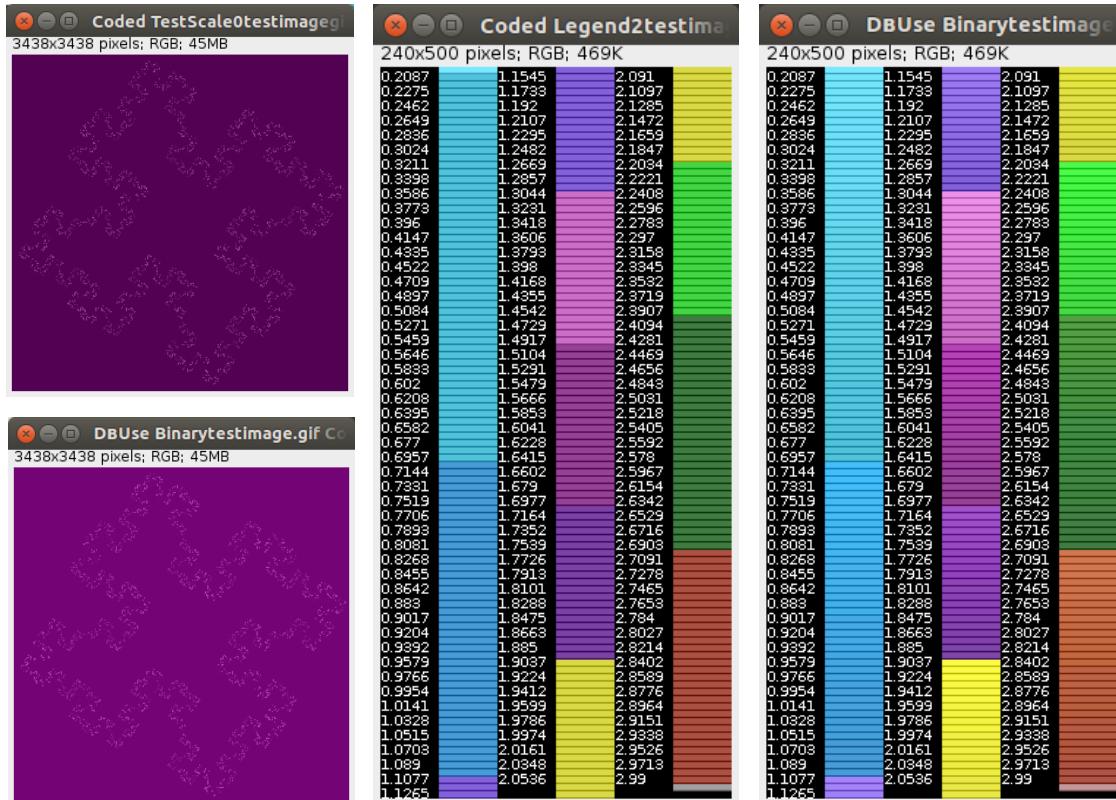


Figure 24: Analyzed Image

## 8. Further Analyze the image according to needs.

## Image Classification

Different types of tumors and diseases are classified separately in this step.

**Recommended tool** *LA-iMageS*

### LA-iMageS [10]

**A software for elemental distribution bioimaging:** Provides easy installation on arbitrary platforms and a simple user interface.

1. Run the following commands in the terminal:

```
$ cd /usr/local/LA-iMageS-1.1.5  
$ ./run.sh
```

2. Select data directory of downloaded tutorial datasets:

*/coin*

*/seed*

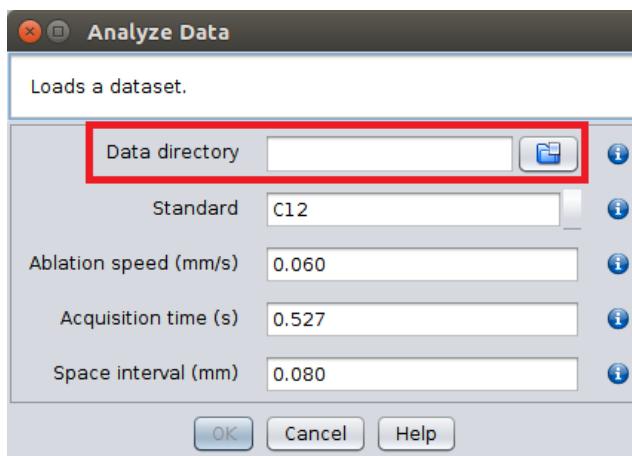


Figure 25: Load Dataset screen

3. Select OK.

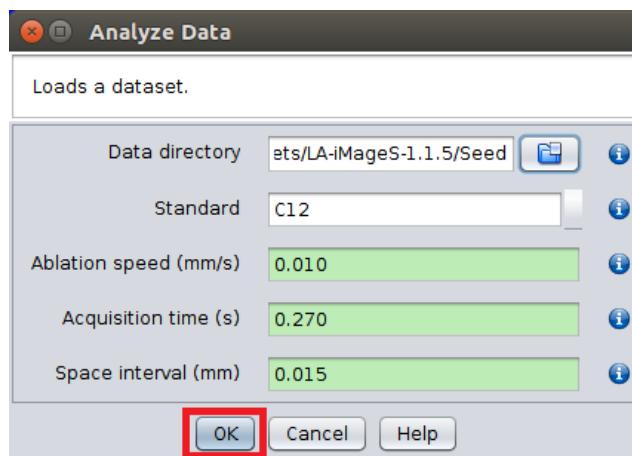


Figure 26: Analyze the Images

**NOTE** Wait for the image to be processed.

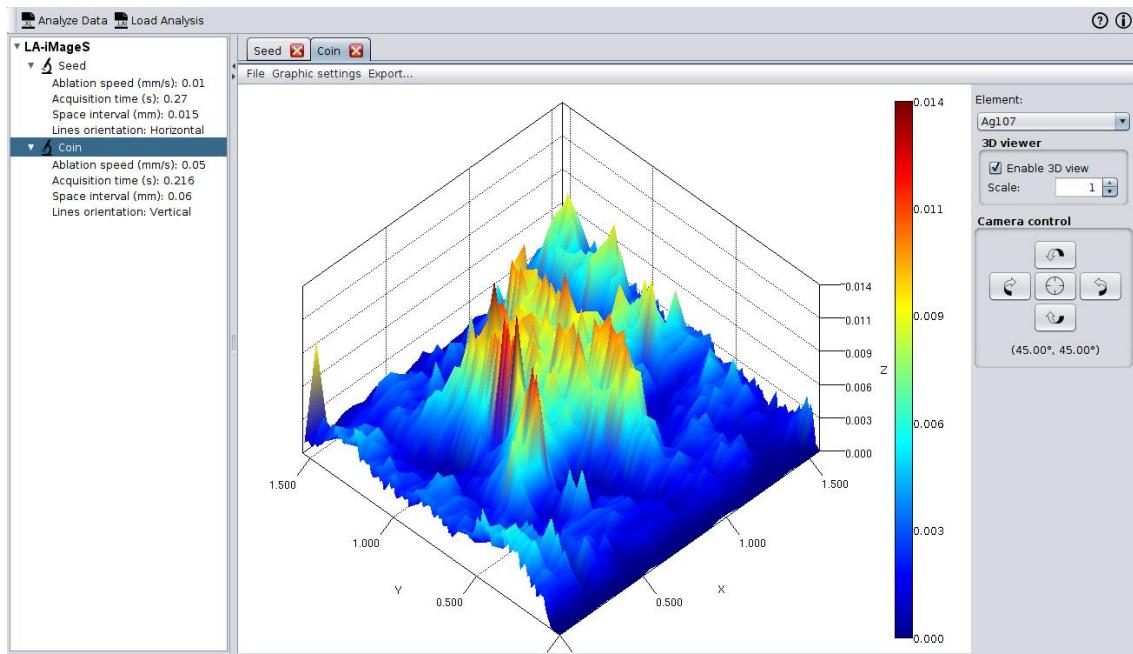


Figure 27: Analysed Image

#### 4. Further Analyze the image.

## Image Mapping

Optimization of images according to T1-weighted and T2-weighted images standards. In T1 images only fat is bright while in T2 images both fat and water is bright. Both of these standards are used for specific analysis.

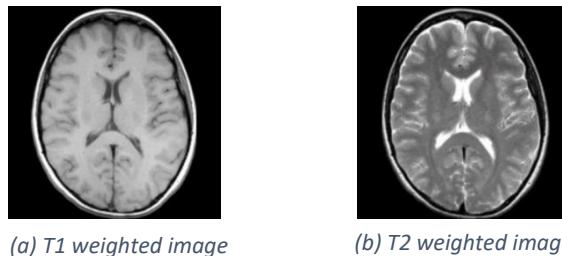


Figure 28: Comparison Between T1 and T2 weighted images

---

### Recommended tool **MRI Processor**

#### MRI Processor [11]

**Computes parametric maps in magnetic resonance (MR) images:** MRI Processor provides different mechanisms to distinguish different tissues and disease processes.

**NOTE** *It is an ImageJ plugin, and it is executed from ImageJ.*

1. Run these commands in the terminal:

```
$ cd ~/Package/Plugins
```

**NOTE** *This is the resource folder of aXonica-master where the tools are downloaded.*

```
$ sudo cp mri_processor_.jar /usr/local/ImageJ/plugins
```

**NOTE** *Provide the password for administrator.*

```
$ cd /usr/local/ImageJ
```

```
$ ./ImageJ
```

2. File → Open

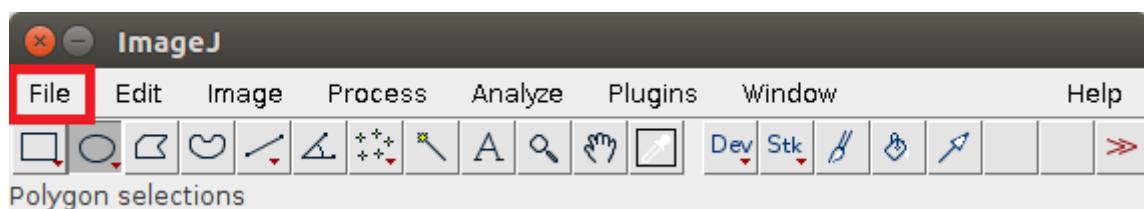


Figure 29: ImageJ Welcome Screen

3. Select the downloaded tutorial datasets: *subject01.nii*

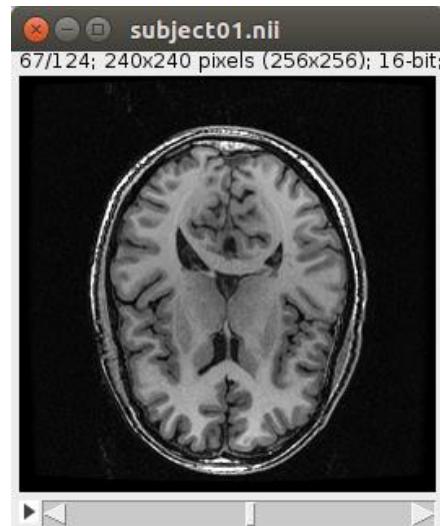


Figure 30: Sample Image

4. Select Plugins → MRI Processor

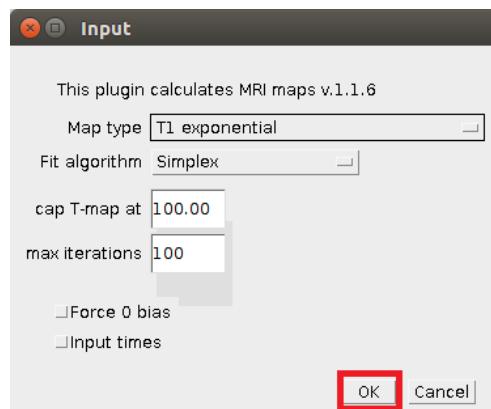


Figure 31: Mapping selection


5. Select Map type as T1 exponential. Click OK

6. Following Mapped Images appear:



Figure 32: Mapped Images

7. Following are the Logs for Image Mapping:

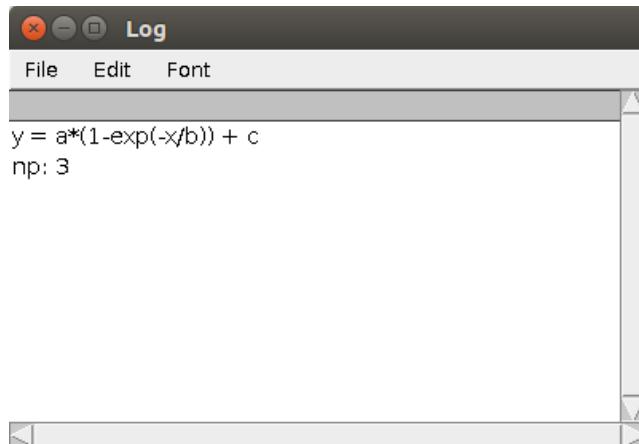

The image shows a window titled "Log" with a dark header bar containing standard window controls (close, minimize, maximize) and the title "Log". Below the header is a menu bar with "File", "Edit", and "Font" options. The main content area is a scrollable text box displaying the following code:

```
y = a*(1-exp(-x/b)) + c  
np: 3
```

The text box has a vertical scrollbar on the right side. At the bottom of the window is a horizontal scrollbar.

Figure 33: Mapping Logs

8. Further analyze the Image

# Data Management and Annotation

---

## In this Chapter

We will learn about:

- Data Management and Annotation
  - Steps for Data Management and Annotation
- 

## Introduction

Refers to the transmission of DICOM image file over networks as well as extracting data from DICOM files and converting it to other formats.

Data Management and Annotation has two significant steps:

- Image Format Management
- DICOM File Management

## Image Format Management

Concerns the transmission of DICOM image file over networks.

**Recommended tool** *MRICron*

### MRICron [12]

**An image viewer for neuroimaging data:** MRICron is a platform able to support multiple layers, draw an identified region of brain injury, view data volume rendering and computes statistical results.

1. Run the following commands in the terminal:

```
$ cd/usr/local/mricron_lx  
$ ./mricron
```

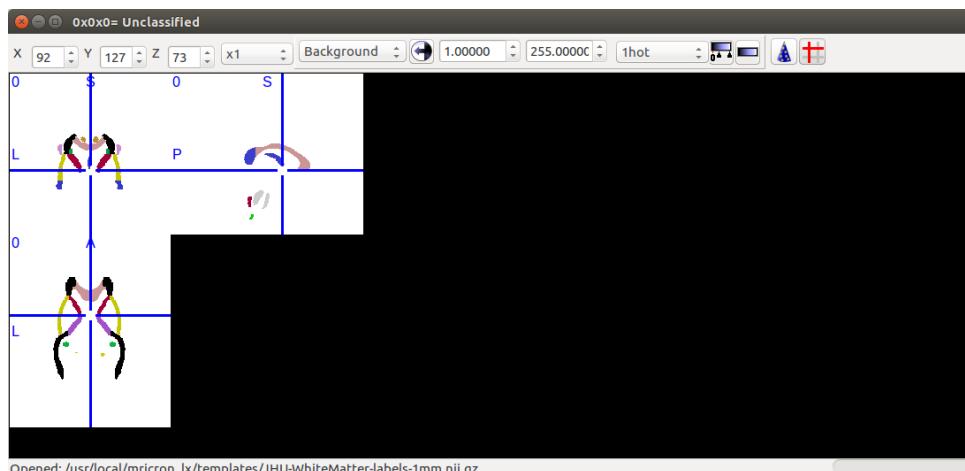


Figure 1: MRICron Welcome screen

2. Further Analyze the image

**NOTE** It is a demo software and it is only used to set the Data Image Formats and test them using a Reconstruction environment.

## DICOM File Management

Deals with the extraction of data from DICOM files and their conversion into other formats.

**Recommended tool WEASIS**

### WEASIS [13]

**Allows users to view clinical images:** WEASIS is a software consisting of a multipurpose web-based viewer and dealing with digital imaging and communications in medicine (DICOM) dataset.

1. Run the following commands in the terminal:

```
$ cd /usr/local/weasis  
$ ./viewer_linux.sh
```

2. File → Open → DICOM



Figure 2: WEASIS Welcome screen

**3. Provide path for the DICOM file.**

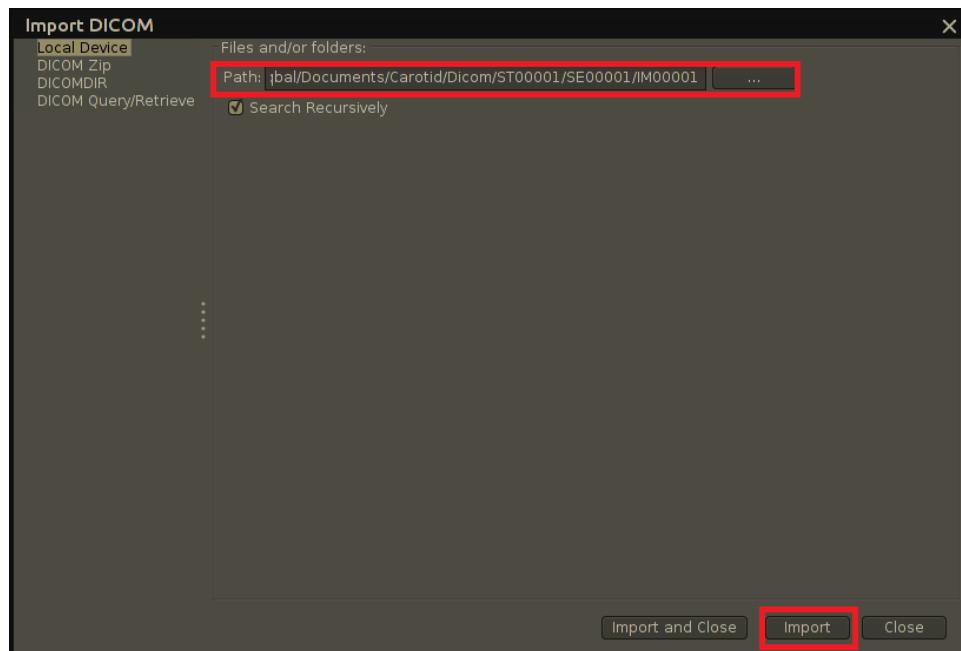


Figure 3: Import DICOM

**4. Select "Import".**

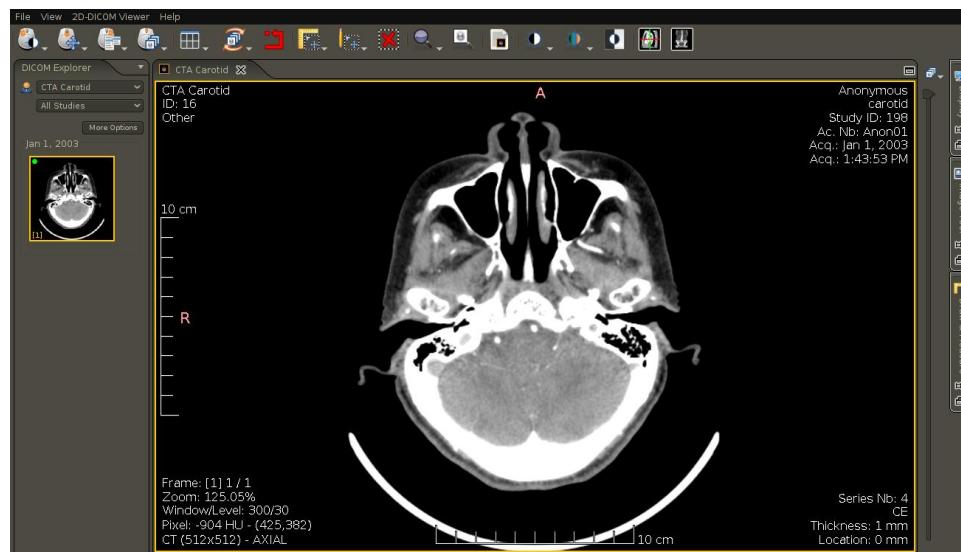


Figure 4: DICOM Image

**5. Further Analyze the image.**

# Post Installation

## Post Installation of aXonica software

After the installation of aXonica is complete you can use the installed tools as follows:

**NOTE** *aXonica\_UI Application is also provided for easier approach..*

### Anatomist

It is installed and used with BrainVISA.

```
$ cd /brainvisa-4.5.0  
$ ./anatomist
```

### Dipy

This tool is a python setting for **MRI** database.

### BoneJ

It is a plugin for ImageJ. It is installed in: Package/Plugins folder. Copy the .jar file and paste in the plugins folder of ImageJ for usage.

```
$ cd /usr/local/ImageJ  
$ ./ImageJ
```

### ImageVIS 3D

```
$ cd /usr/local/ImageVis3D-3.1.0/  
$ ./ImageVis3D
```

### AMIDE

```
$ export UBUNTU_MENUPROXY=0  
$ amide
```

### PID

This is a command line interface tool. Its user manual is provided.

## TORTOISE

```
$ cd /usr/local/TORTOISE_V3.1.0/TORTOISE_GUI  
$ gksu python TortoiseGui.py
```

## normalizeFOV

This is a command line interface tool. Its user manual is provided.

## FastFilter3D

It is a plugin for ImageJ. It is installed in: Package/Plugins folder. Copy the .jar file and paste in the plugins folder of ImageJ for usage.

```
$ cd /usr/local/ImageJ  
$ ./ImageJ
```

## ITK-Snap

```
$ cd /usr/local/itksnap-3.6.0-20170401-Linux-x86_64  
$ cd bin/  
$ ./itksnap
```

## HeteroscedasticfMRI

This is a command line interface tool. Its user manual is provided.

## pyClusterROI

This is a command line interface tool. Its user manual is provided.

## BrainSeg3D

```
$ cd /usr/local/brainseg3d  
$ ./BrainSeg3D
```

## CMP-BIA

It is a plugin for ImageJ. It is installed in: Package/Plugins folder. Copy the .jar file and paste in the plugins folder of ImageJ for usage

```
$ cd /usr/local/ImageJ  
$ ./ImageJ
```

## ITK

This is a command line interface tool. Its user manual is provided.

## E-Snake

It is a plugin for ImageJ. It is installed in: Package/Plugins folder. Copy the .jar file and paste in the plugins folder of ImageJ for usage

```
$ cd /usr/local/ImageJ  
$ ./ImageJ
```

## mincBEAST

This is command line interface tool. Its user manual is provided.

## Segmentator

This is command line interface tool. Its user manual is provided.

## Mango

```
$ cd /usr/local/Mango  
$ ./mango
```

## bUnwarpJ

It is a plugin for ImageJ. It is installed in: Package/Plugins folder. Copy the .jar file and paste in the plugins folder of ImageJ for usage

```
$ cd /usr/local/ImageJ  
$ ./ImageJ
```

## BrainVISA

```
$ cd /brainvisa-4.5.0  
$ ./anatomist
```

## The Virtual Brain

```
$ cd /usr/local/TVB_Distribution/bin/  
$ ./tvb_start.sh
```

## CBS Tools

It is a plugin for LA-iMageS.

```
> cd /usr/local/LA-iMageS-1.1.5  
> chmod +x run.sh  
> ./run.sh
```

## HegaEtAI2017

This is a command line interface tool. Its user manual is provided.

## NEURON

This is a command line interface tool. Its user manual is provided.

## SliceMap

It is a plugin for ImageJ. It is installed in: Package/Plugins folder. Copy the .jar file and paste in the plugins folder of ImageJ for usage

```
$ cd /usr/local/ImageJ  
$ ./ImageJ
```

## Nengo

This is command line interface tool. Its user manual is provided.

## neuroConstruct

```
$ cd /usr/local/neuroConstruct_1.6.0  
$ gksu bash neuroConstruct_1.6.0
```

## MedInria

```
$ cd /usr/local/medinria-2.2.3-Linux-x86_64  
$ cd bin/  
$ ./medInria_launcher.sh
```

## Brian

This is command line interface tool. Its user manual is provided.

## Time Domain decoding

This is command line interface tool. Its user manual is provided.

## Fraclac

It is a plugin for ImageJ. It is installed in: Package/Plugins folder. Copy the .jar file and paste in the plugins folder of ImageJ for usage

```
$ cd /usr/local/ImageJ  
$ ./ImageJ
```

## Gwyddion

```
$ gwyddion
```

## LA-iMageS

```
$ cd /usr/local/LA-iMageS-1.1.5  
$ chmod +x run.sh  
$ ./run.sh
```

## Brainstorm

```
$ cd /usr/local/brainstorm3/bin/R2015b/  
$ gksu bash brainstorm3.command
```

## ImageJ

```
$ cd /usr/local/ImageJ  
$ ./ImageJ
```

## Nipype

This tool is python setting for MRI database.

## MRI Processor

It is a plugin for ImageJ. It is installed in: Package/Plugins folder. Copy the .jar file and paste in the plugins folder of ImageJ for usage

```
$ cd /usr/local/ImageJ  
$ ./ImageJ
```

## NiBabel

This tool is python setting for MRI database.

## MRICron

```
$ cd /usr/local/mricron_lx  
$ ./mricron
```

## MIPAV

```
$ cd /usr/local/MIPAV  
$ ./mipav
```

## WEASIS

```
$ cd /usr/local/weasis  
$ ./viewer-linux.sh
```

# Software Uninstallation Process

## Uninstalling aXonica software

1. Uninstaller is downloaded with the Installer of aXonica.
2. Run the following commands on terminal:  
`> chmod +x aXonica_uninstall_enUS`  
`> sudo bash aXonica_uninstall_enUS`
3. Uninstallation wizard of aXonica will start.

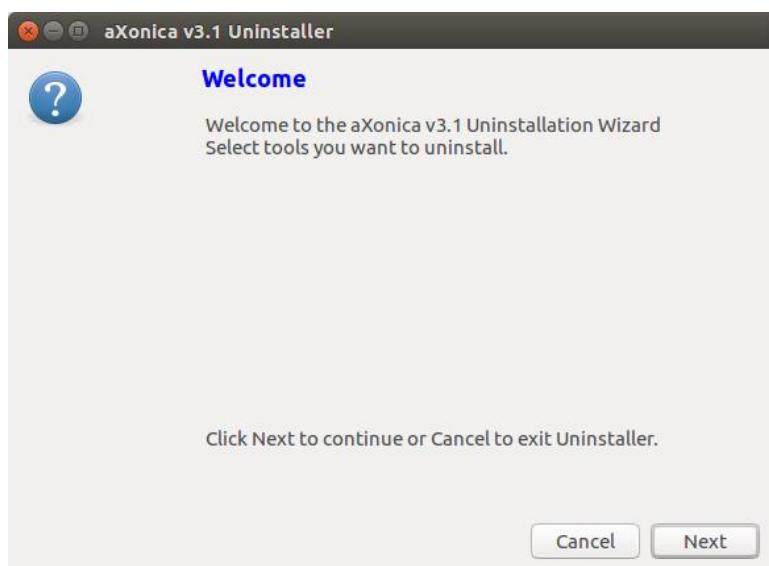


Figure 1: Uninstallation Welcome screen

**NOTE** Software which are not installed also appear in the uninstaller selection screens. Selecting/Deselecting them has no effect on uninstallation process.

4. Selection screen for Preprocessing tools is displayed. Select the required tools and Click **Next**.

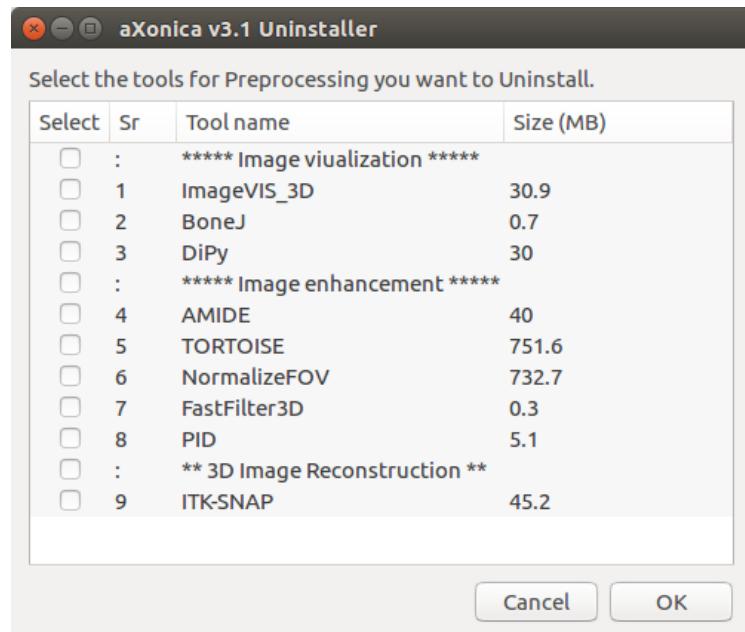


Figure 2: Preprocessing tools selection screen

5. Selection screen for Processing tools is displayed. Select the required tools and Click **Next**.

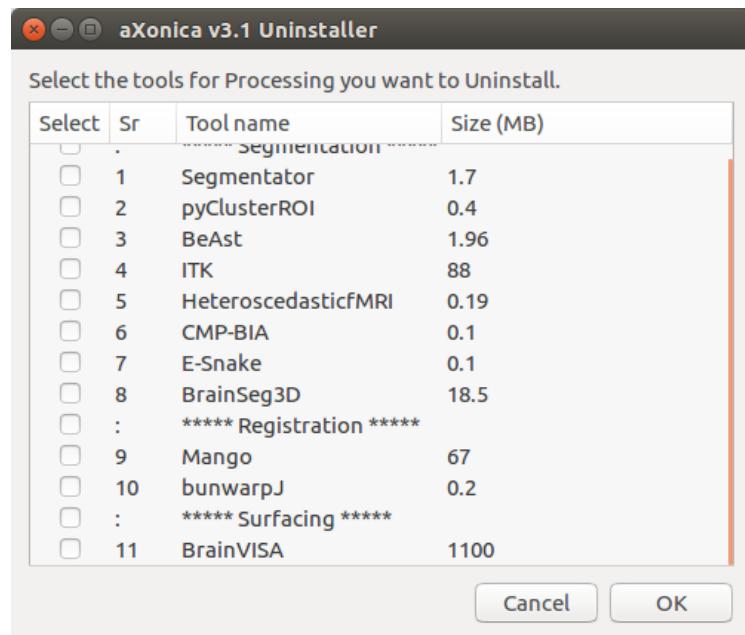


Figure 3: Processing tools selection screen

6. Selection screen for Structural Analysis tools is displayed. Select the required tools and Click **Next**.

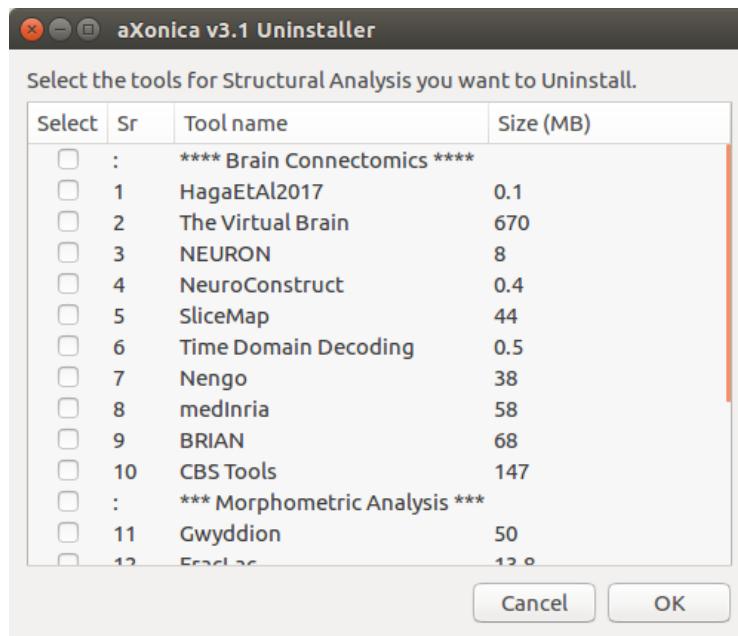


Figure 4: Structural Analysis tools selection screen

7. Selection screen for Data Management and Annotation tools is displayed. Select the required tools and Click **Next**.

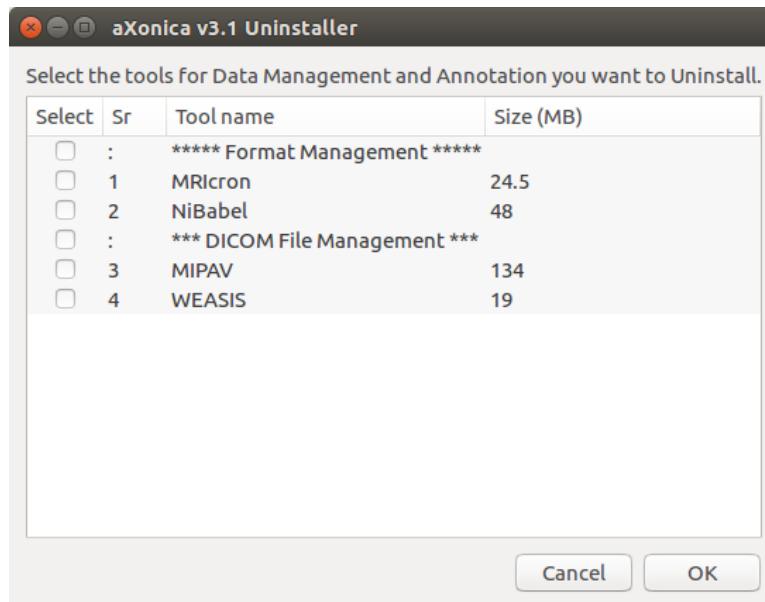


Figure 5: Data Management and Annotation selection screen

8. Click **Next** to proceed with the uninstallation.

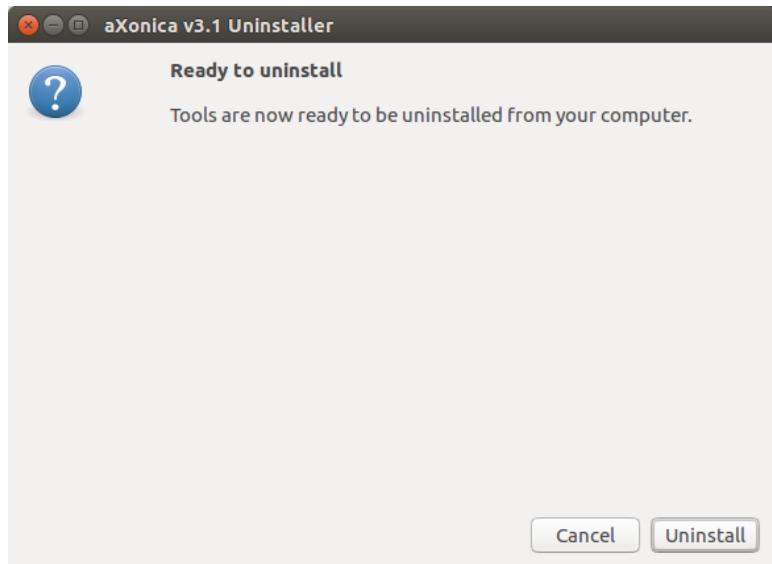


Figure 6: Ready to uninstall screen

23. Tools uninstallation will **continue**.

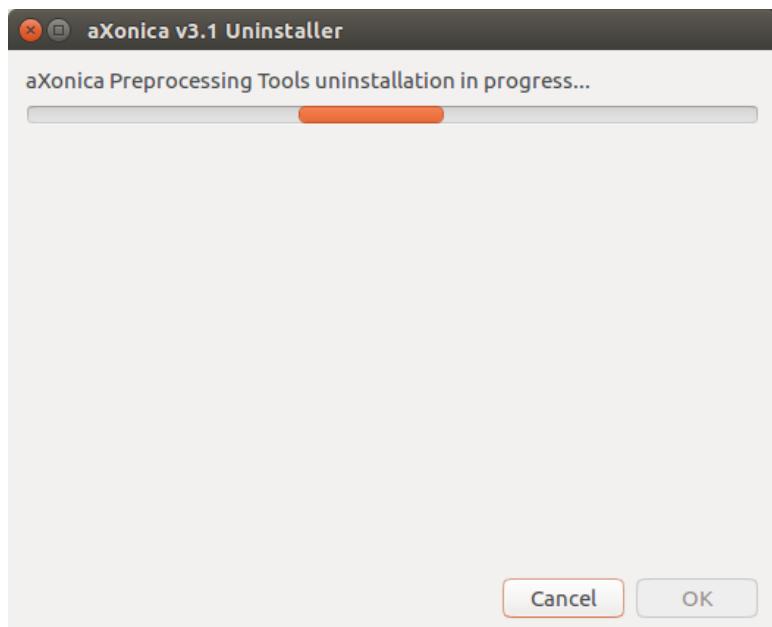


Figure 7: Tools Uninstallation progress screen

9. Uninstallation of aXonica is now finished. Click **Finish** to continue.

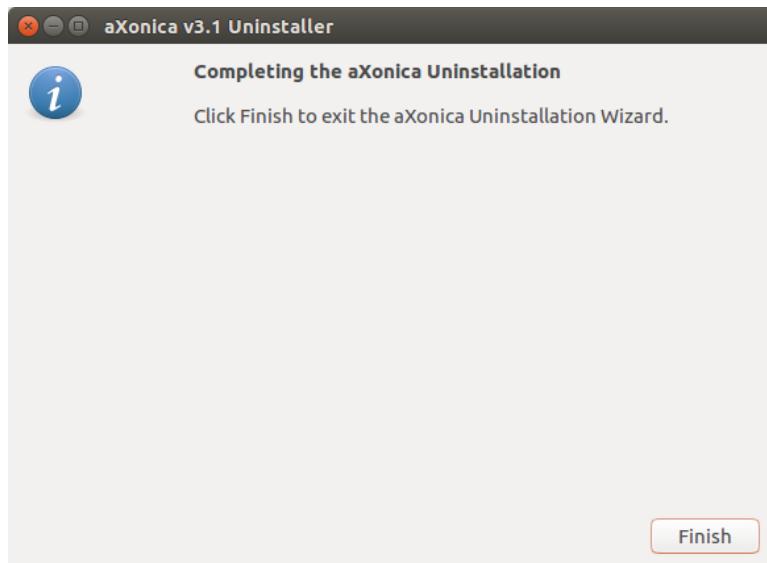


Figure 8: Finalizing Uninstallation screen

# References

- [1] T. Fogal and J. H. Kruger, "Tuvok, an Architecture for Large Scale Volume Rendering," in *VMV*, 2010.
- [2] A. M. Loening and S. S. Gambhir, "AMIDE: a free software tool for multimodality medical image analysis," *Molecular imaging*, vol. 2, p. 15353500200303133, 2003.
- [3] P. A. Yushkevich, J. Piven, H. C. Hazlett, R. G. Smith, S. Ho, J. C. Gee and G. Gerig, "User-guided 3D active contour segmentation of anatomical structures: significantly improved efficiency and reliability," *Neuroimage*, vol. 31, pp. 1116-1128, 2006.
- [4] Z. Lesjak, A. Galimzianova, B. Likar, F. Pernus and Z. Spiclin, "Increased accuracy and reproducibility of MS lesion volume quantification by using publicly available BrainSeg3D image analysis software," in *MULTIPLE SCLEROSIS JOURNAL*, 2015.
- [5] J. L. Lancaster, A. R. Laird, S. B. Eickhoff, M. J. Martinez, P. M. Fox and P. T. Fox, "Automated regional behavioral analysis for human brain images," *Frontiers in neuroinformatics*, vol. 6, p. 23, 2012.
- [6] Y. Cointepas, J.-F. Mangin, L. Garnero, J.-B. Poline and H. Benali, "BrainVISA: software platform for visualization and analysis of multi-modality brain data," *Neuroimage*, vol. 13, p. 98, 2001.
- [7] N. Toussaint, J.-C. Souplet, P. Fillard and others, "MedINRIA: medical image navigation and research tool by INRIA," in *Proc. of MICCAI*, 2007.
- [8] J. Schindelin, C. T. Rueden, M. C. Hiner and K. W. Eliceiri, "The ImageJ ecosystem: An open platform for biomedical image analysis," *Molecular reproduction and development*, vol. 82, pp. 518-529, 2015.
- [9] A. Karperien, "FracLac for ImageJ," *Charles Sturt University*, 2013.
- [10] H. Lopez-Fernandez, G. S. Pessoa, M. A. Z. Arruda, J. L. Capelo-Martinez, F. Fdez-Riverola, D. Glez-Peña and M. Reboiro-Jato, "LA-iMageS: a software for elemental distribution bioimaging using LA--ICP--MS data," *Journal of cheminformatics*, vol. 8, p. 65, 2016.
- [11] D. Prodanov and K. Verstreken, "Automated segmentation and morphometry of cell and tissue structures. Selected algorithms in imageJ," in *Molecular Imaging*, InTech, 2012.
- [12] S. R. Kesler, J. S. Kent and R. O'Hara, "Prefrontal cortex and executive function impairments in primary breast cancer," *Archives of neurology*, vol. 68, pp. 1447-1453, 2011.
- [13] G. Valeri, F. A. Mazza, S. Maggi, D. Aramini, L. La Riccia, G. Mazzoni and A. Giovagnoni, "Open source software in a practical approach for post processing of radiologic images," *La radiologia medica*, vol. 120, pp. 309-323, 2015.